

# **XPS data analysis code: Ctrl+Q**

*- efficient and effective analysis ever experienced -*

Hideki NAKAJIMA

*heitler@gmail.com*


*github: heitler/xps-excel-macro*

23 Aug. 2016


# Brief

## Plot with binding energy


### Identify elements and states

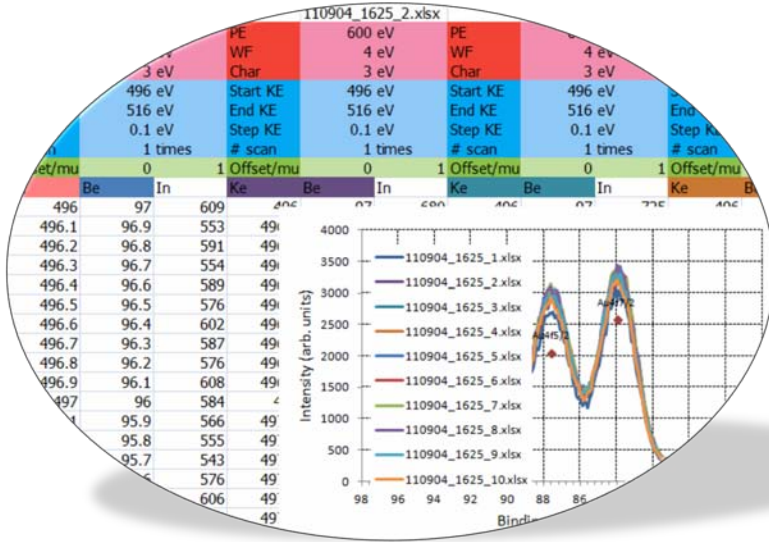

# Energy and intensity calibrations

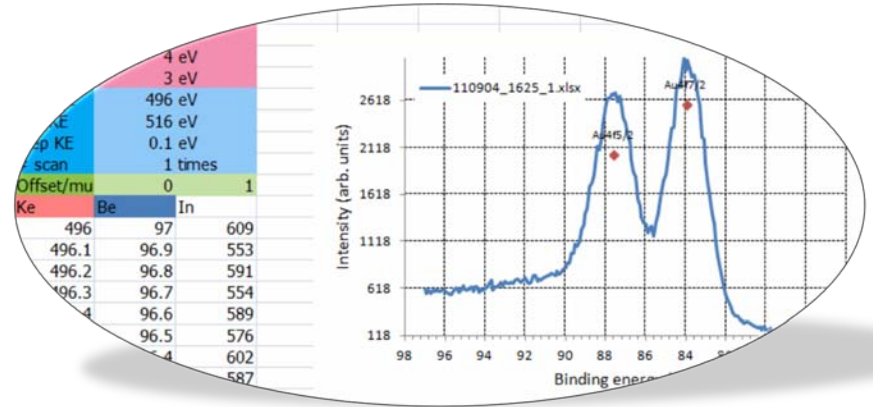

## Curve fitting and BG subtraction

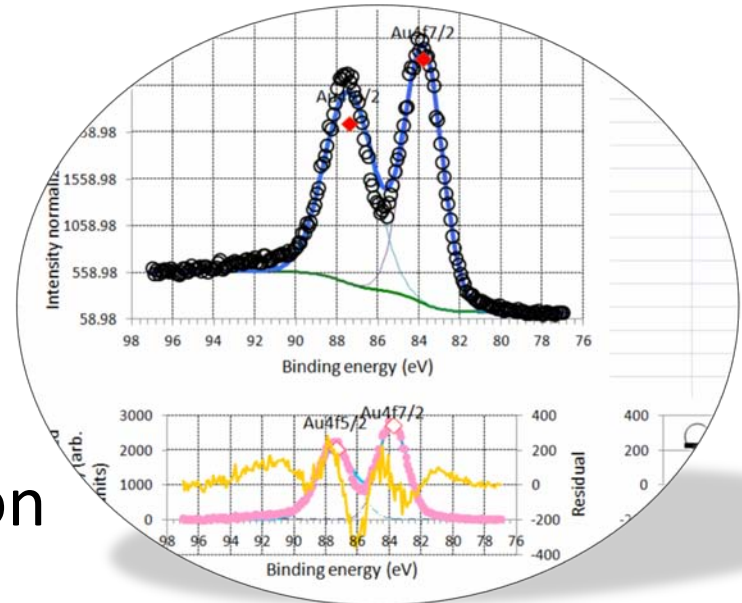



# Preface


- **Ctrl+Q** is the first comprehensive data-analysis code for **SR-XPS** and **XAS** working on Windows Office **Excel 2007 or later**.
- Ctrl+Q will analyze the data by just pressing the shortcut key: **Ctrl + q**.
- It works great for the number of publications.
  - Carbon, JPAP, ASS, ACS AMI, SAB etc. at BL3.2Ua in SPL
- **GitHub: xps-excel-macro** for updated info.

 [heitler / xps-excel-macro](#)

 Watch  Star  Fork

 Code

 Issues 0

 Pull requests 0

 Wiki

 Pulse

 Graphs

 Settings


The VBA codes to be pasted in the personal macro workbook (Personal.xlsb) in Windows Office Excel 2007 - 2016 Visual Basic Editor (VBE) work with Solver for spectral analysis and curve fitting of the synchrotron-based soft x-ray photoemission (XPS) and absorption (XAS) spectra.

— Edit

# Data preparation

	A	B
1	KE/eV	Scan#1
2	496	19.69872
3	496.1	17.95129
4	496.2	19.18263
5	496.3	18.03784
6	496.4	19.03867
7	496.5	18.80417
8	496.6	18.22572


- Two-column data in the Tab-delimited format texts or spreadsheets are analyzed in the code.
- Sheet is named after the workbook name.
- Syntax in A1 cell has to be chosen as follows.

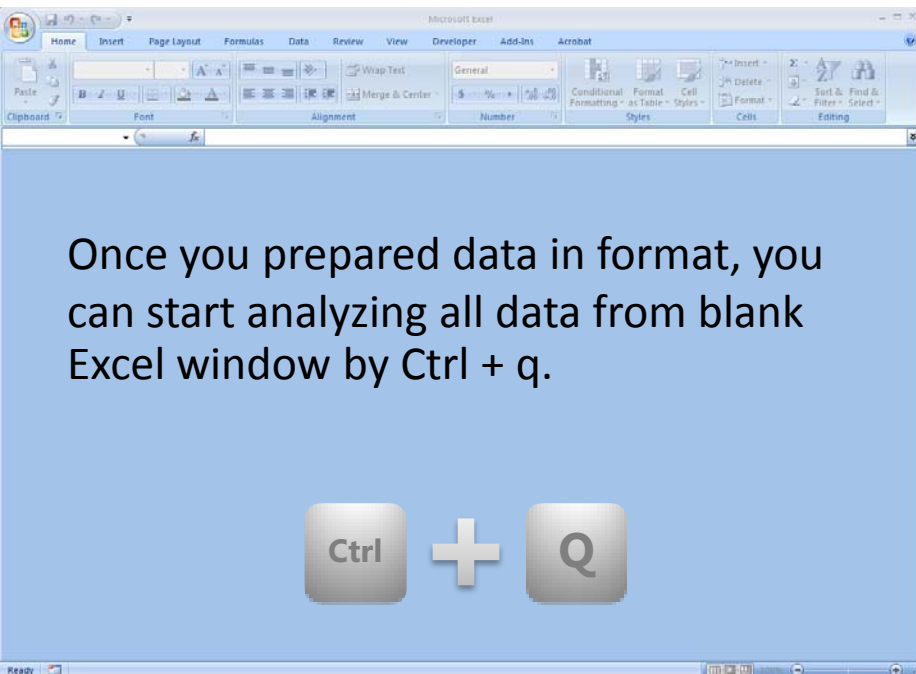

Syntax in A1 cell	Data in column A	Data in column B
KE/eV	KE	XPS
BE/eV	BE	XPS
PE/eV	PE	XAS
GE/eV	PE	G scan
AE/eV	EE	Auger
QE/eV	mass	Q-mass
ME/eV	Position	Any

# Installation of the code in Excel

- Generate **Personal Workbook Macro**
- Install **Solver** Add-in and setup its DLL
- Setup the **Shortcut key** in Macro Option

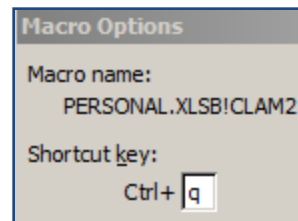

- Compile your own **database** for BE and sensitivities on UD.xlsx


XPS sheet

	A	B	C	D
1	Element	Orbit	BE(eV)	ASF
2	C	1s	284.6	1
3	O	1s	532	2.93

AES sheet

	A	B	C	D
1	Element	Auger	KE(eV)	RSF
2	C	KLL	266	0.6
3	O	KLL	506	0.96


# C1s sample csv

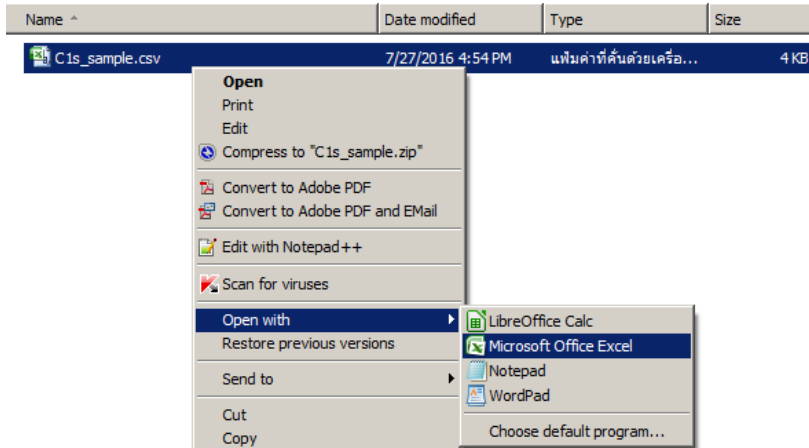


	A	B
1	KE/eV	PE: 600 eV
2	290	1709
3	290.1	1798
4	290.2	1845
5	290.3	1722
6	290.4	1729
7	290.5	1699
8	290.6	1727
9	290.7	1693
10	290.8	1647
11	290.9	1684
12	291	1650



600 eV of photon energy

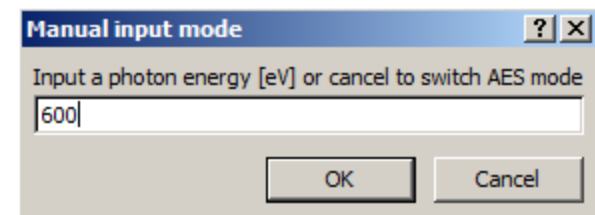

XPS: Carbon 1s spectrum

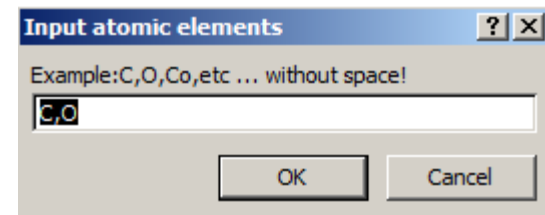

# Sample file with multiple scans

File Explorer context menu for **Au4f\_sample.csv**:

- Open
- Print
- Edit
- Compress to "Au4f\_sample.zip"
- Convert to Adobe PDF
- Convert to Adobe PDF and Email
- Edit with Notepad ++
- Scan for viruses
- Open with
  - LibreOffice Calc
  - Microsoft Office Excel**
  - Notepad
  - WordPad
- Restore previous versions
- Send to
- Cut
- Copy
- Create shortcut
- Delete
- Rename
- Properties

Microsoft Excel spreadsheet data (approximate):

	A	B	C	D	E	F	G	H
1	KE/eV	110904_16	KE/eV	110904_1	KE/eV	110904_16	KE/eV	110904_1
2	496	609	496	680	496	725	496	
3	496.1	553	496.1	660	496.1	680	496.1	
4	496.2	591	496.2	653	496.2	660	496.2	
5	496.3	554	496.3	645	496.3	660	496.3	
6	496.4	589	496.4	659	496.4	693	496.4	
7	496.5	576	496.5	687	496.5	660	496.5	

Shortcut keys: **Ctrl + Q**

File Explorer list:

Name	Date	Type
110904_1625_1.txt	6/17/2016 1:29 PM	TXT File
110904_1625_2.txt	6/17/2016 1:29 PM	TXT File
110904_1625_3.txt	6/17/2016 1:29 PM	TXT File
110904_1625_4.txt	6/17/2016 1:29 PM	TXT File
110904_1625_5.txt	6/17/2016 1:29 PM	TXT File
110904_1625_6.txt	6/17/2016 1:29 PM	TXT File
110904_1625_7.txt	6/17/2016 1:29 PM	TXT File
110904_1625_8.txt	6/17/2016 1:29 PM	TXT File
110904_1625_9.txt	6/17/2016 1:29 PM	TXT File
110904_1625_10.txt	6/17/2016 1:29 PM	TXT File
110904_1625_11.txt	6/17/2016 1:29 PM	TXT File

Microsoft Excel dialog box: Data were exported in the text files. [OK]

- Download sample.txt including 100 spectra.
- Open with Excel as shown above figure.
- Export 100 text files by shortcut keys Ctrl+Q.

# Sample file with single scan

Name ^	Date modified	Type	Size
110904_1625_1.txt	7/27/2016 5:12 PM	TXT File	3 KB
110904_1625_2.txt		TXT File	3 KB
110904_1625_3.txt		TXT File	3 KB
110904_1625_4.txt		TXT File	3 KB
110904_1625_5.txt		TXT File	3 KB
110904_1625_6.txt		TXT File	3 KB
110904_1625_7.txt		TXT File	3 KB
110904_1625_8.txt		TXT File	3 KB
110904_1625_9.txt		TXT File	3 KB
110904_1625_10.txt		TXT File	3 KB
110904_1625_11.txt		TXT File	3 KB
110904_1625_12.txt		TXT File	3 KB
110904_1625_13.txt		TXT File	3 KB
110904_1625_14.txt		TXT File	3 KB
110904_1625_15.txt		TXT File	3 KB
110904_1625_16.txt		TXT File	3 KB
110904_1625_17.txt		TXT File	3 KB
110904_1625_18.txt		TXT File	3 KB

Open

Edit

Compress to "110904\_1625\_1.zip"

Convert to Adobe PDF

Combine supported files in Acrobat...

Edit with Notepad++

Scan for viruses

Open with

Restore previous versions

Send to

Cut

Copy

Create shortcut

Delete

Rename

Properties

Adobe Illustrator CS6

LibreOffice Calc

LibreOffice Writer

Mery

Microsoft Office Excel

Notepad

Notepad++ : a free (GNU) source co

WordPad

Choose default program...



	A	B
1	KE/eV	110904_16
2	496	609
3	496.1	553
4	496.2	591
5	496.3	554
6	496.4	589
7	496.5	576
8	496.6	602
9	496.7	587
10	496.8	576
11	496.9	608
12	497	594



# Single data to plot in Graph sheet

	A	B
1	KE/eV	110904.16
2	496	609
3	496.1	553
4	496.2	591
5	496.3	554
6	496.4	589
7	496.5	576
8	496.6	602
9	496.7	587
10	496.8	576
11	496.9	608
12	497	584

KE/eV represents a kinetic energy.



600 eV of photon energy



**Manual input mode** ? X

Input a photon energy [eV] or cancel to switch AES mode

OK Cancel

Sample: Gold metal foil



**Input atomic elements** ? X

Example: C,O,Co,etc ... without space!

OK Cancel

Note that Au4f sensitivity is not in database, so you have to input factors in database in a way below.

Au 4f BE and ASF  
in XPS sheet on UD.xlsx



	A	B	C	D
1	Element	Orbit	BE(eV)	ASF
2	C	1s	284.6	1
3	O	1s	532	2.93
4	Au	4f5/2	87.6	7.54
5	Au	4f7/2	84	9.58


[illegible]

A plot of Intensity (arb. units) versus Kinetic energy (eV). The x-axis ranges from 496 to 516 eV, and the y-axis ranges from 118 to 2618. A red line shows the experimental data, which exhibits a broad peak around 508 eV and a sharper peak around 512 eV. Two specific peaks are labeled: Au 4f5/2 at approximately 508.5 eV and Au 4f7/2 at approximately 512.5 eV. An orange box with the text "KE plot" is overlaid on the right side of the plot.

# Reference spectrum

Grating	0 lines/mm		
PE	600 eV		
WF	4 eV		
Char	3 eV		
Start KE	496 eV		
End KE	516 eV		
Step KE	0.1 eV		
# scan	1 times		
Offset/mu	0	1	
Ke	Be	In	
496	97	609	
496.1	96.9	553	
496.2	96.8	591	
496.3	96.7	554	
496.4	96.6	589	
496.5	96.5	576	
496.6	96.4	602	
496.7	96.3	587	
496.8	96.2	576	

Standard peak BE should be well-known and assumed to be identical with exp.

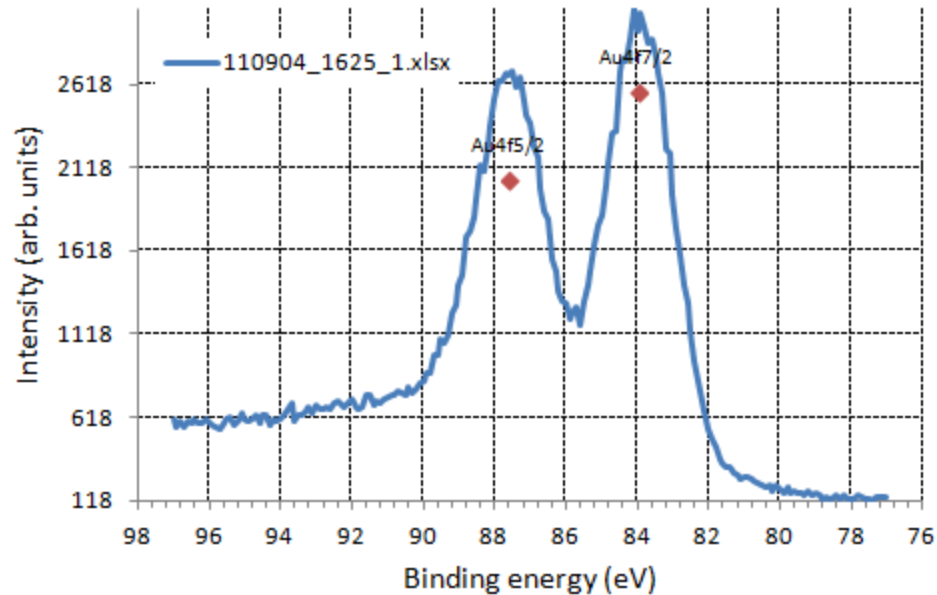


Next, apply these factors to another by “debug” in A1 cell then



	A	B	C
1	debug	0 lines/mm	
2	PE	600 eV	
3	WF	4 eV	
4	Char	3 eV	
5	Start KE	496 eV	
6	End KE	516 eV	
7	Step KE	0.1 eV	
8	# scan	1 times	
9	Offset/mu	0	1
10	Ke	Be	In

# Multiple-file selection

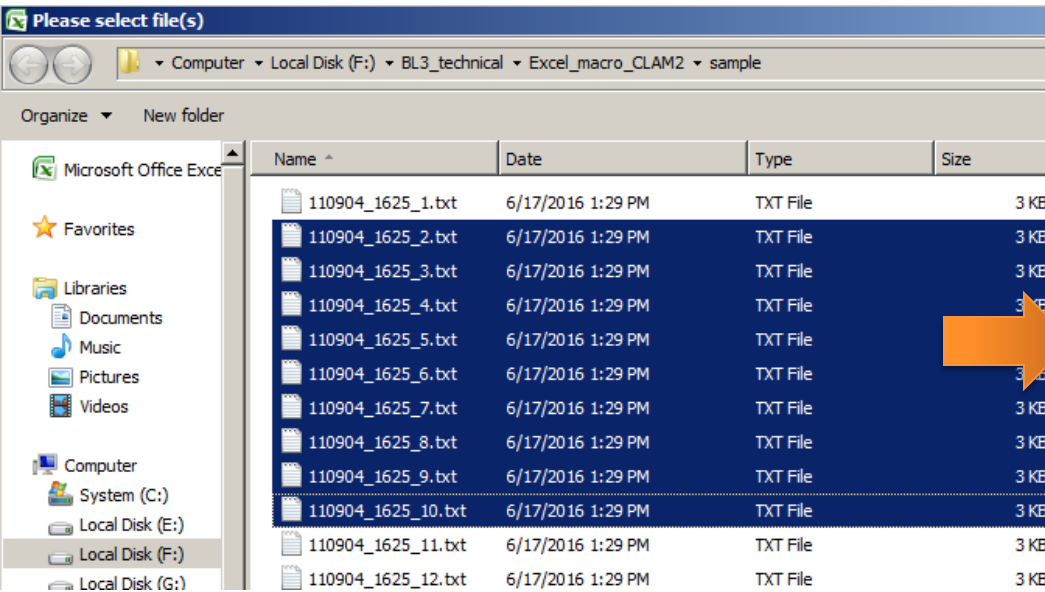

Name	Date	Type	Size
110904_1625_1.txt	6/17/2016 1:29 PM	TXT File	3 KB
110904_1625_1.xlsx	6/17/2016 1:55 PM	แผ่นงาน Microsoft O...	53 KB
110904_1625_2.txt	6/17/2016 1:29 PM	TXT File	3 KB
110904_1625_2.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	52 KB
110904_1625_3.txt	6/17/2016 1:29 PM	TXT File	3 KB
110904_1625_3.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	12 KB
110904_1625_4.txt	6/17/2016 1:29 PM	TXT File	3 KB
110904_1625_4.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	53 KB
110904_1625_5.txt	6/17/2016 1:29 PM	TXT File	3 KB
110904_1625_5.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	53 KB
110904_1625_6.txt	6/17/2016 1:29 PM	TXT File	3 KB
110904_1625_6.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	52 KB
110904_1625_7.txt	6/17/2016 1:29 PM	TXT File	3 KB
110904_1625_7.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	52 KB
110904_1625_8.txt	6/17/2016 1:29 PM	TXT File	3 KB

Next, compare the spectra by "comp" in D1

	A	B	C	D
1	Grating	0 lines/mm	comp	
2	PE	600 eV		
3	WF	4 eV		
4	Char	3 eV		
5	Start KE	496 eV		
6	End KE	516 eV		
7	Step KE	0.1 eV		
8	# scan	1 times		
9	Offset/mu	0	1	
10	Ke	Be	In	
11	496	97	609	
12	496.1	96.9	553	
13	496.2	96.8	591	
14	496.3	96.7	554	
15	496.4	96.6	589	
16	496.5	96.5	576	
17	496.6	96.4	602	
18	496.7	96.3	587	




Name	Date modified	Type	Size
110904_1625_1.xlsx	6/17/2016 2:06 PM	แผ่นงาน Microsoft O...	177 KB
110904_1625_2.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	53 KB
110904_1625_3.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	53 KB
110904_1625_4.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	53 KB
110904_1625_5.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	53 KB
110904_1625_6.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	52 KB
110904_1625_7.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	52 KB
110904_1625_8.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	52 KB
110904_1625_9.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	52 KB
110904_1625_10.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	53 KB
sample.xlsx	6/17/2016 1:55 PM	แผ่นงาน Microsoft O...	181 KB

# Compared spectra

Grating	0 lines/mm		110904_1625_2.xlsx			110904_1625_3.xlsx			110904_1625_4.xlsx
PE	600 eV		PE	600 eV		PE	600 eV		PE
WF	4 eV		WF	4 eV		WF	4 eV		WF
Char	3 eV		Char	3 eV		Char	3 eV		Char
Start KE	496 eV		Start KE	496 eV		Start KE	496 eV		Start KE
End KE	516 eV		End KE	516 eV		End KE	516 eV		End KE
Step KE	0.1 eV		Step KE	0.1 eV		Step KE	0.1 eV		Step KE
# scan	1 times		# scan	1 times		# scan	1 times		# scan
Offset/mu	0		Offset/mu	0		Offset/mu	0		Offset/mu
Ke	Be	In	Ke	Be	In	Ke	Be	In	Ke

496	97	609
496.1	96.9	553
496.2	96.8	591
496.3	96.7	554
496.4	96.6	589
496.5	96.5	576
496.6	96.4	602
496.7	96.3	587
496.8	96.2	576
496.9	96.1	608
497	96	584
497.1	95.9	566
497.2	95.8	555
497.3	95.7	543
497.4	95.6	576
497.5	95.5	606
497.6	95.4	614
497.7	95.3	575

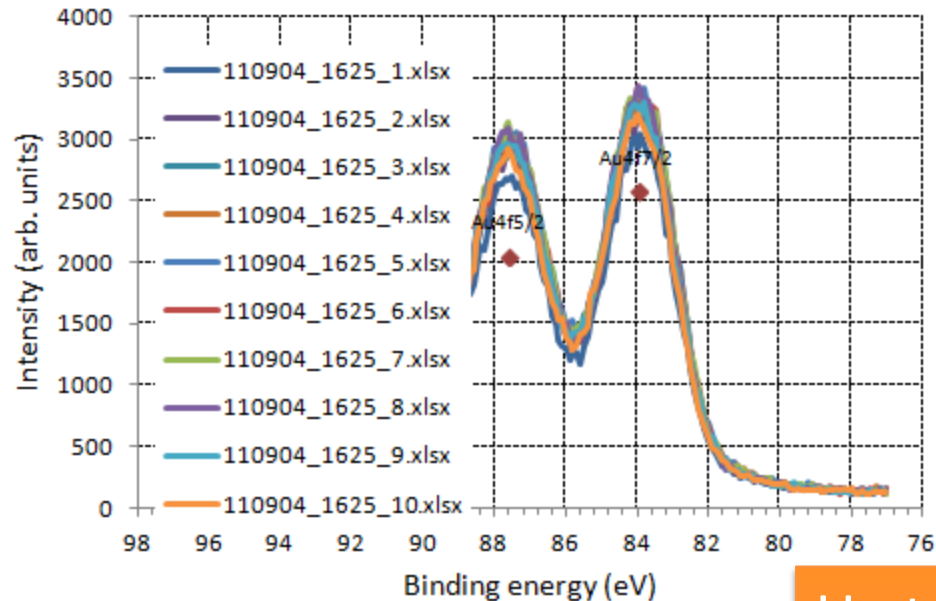

Up to 100 spectra  
for BE&KE plots

# Analysis on spectra compared


- Automatically scaled spectra for comparison
  - “auto” in A1 cell to scale at the both ends
  - “auto[x1:x2,x3:x4]” in A1 cell to scale in the ranges of (x1, x2) for offset (0) and (x3, x4) for multiple (1).
- Annotate legends in each plot
  - “leg” in A1 cell to generate Sample sheet
  - Describe sample or spectrum name in the sheet
  - “leg” in A1 cell to run the code again.
- Normalize spectra with ref. spectrum by “norm” in A1 cell.
  - Second set of data will be selected for ref. and third set to be normalized data.
  - Norm sheet appears from the data normalized.

Every single step needs




	A	B	C		A	B	C	D	E
1	leg	0	lines/mm	1	No.	Name	Sep.	File name	
2	PE	600 eV		2	1	Au1	:	110904_1625_1.xlsx	
3	WF	4 eV		3	2	Au2	:	110904_1625_2.xlsx	
4	Char	3 eV		4	3	Au3	:	110904_1625_3.xlsx	
5	Start KE	496		5	4	Au4	:	110904_1625_4.xlsx	
6	End KE	516		6	5	Au5	:	110904_1625_5.xlsx	
7	Step KE	0.1 eV		7	6	Au6	:	110904_1625_6.xlsx	
8	# scan	1 times		8	7	Au7	:	110904_1625_7.xlsx	
9	Offset/mu	0	1	9	8	Au8	:	110904_1625_8.xlsx	
10	Ke	Be	In	10	9	Au9	:	110904_1625_9.xlsx	
				11	10	Au10	:	110904_1625_10.xlsx	


Sample sheet

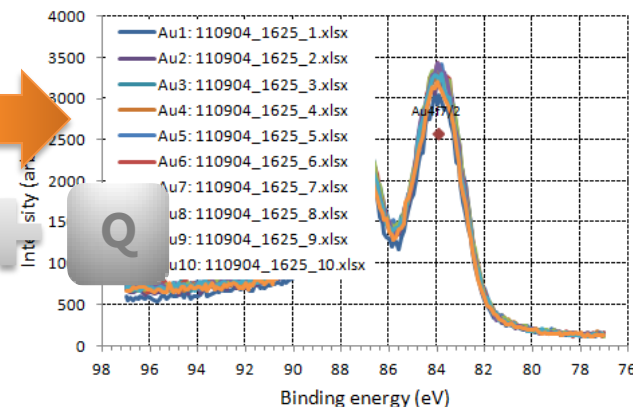

# Fit sheet

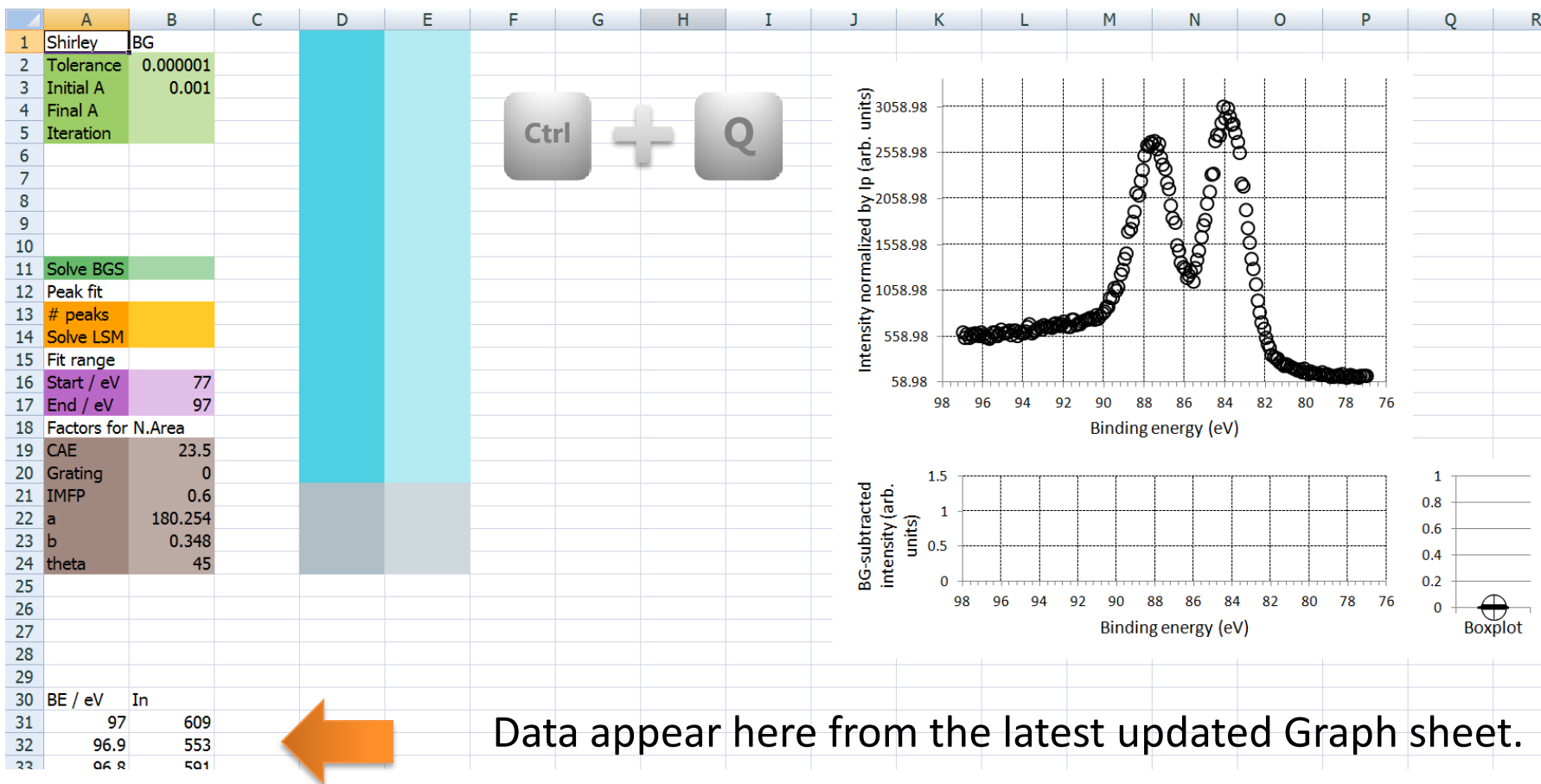


Switch the sheet

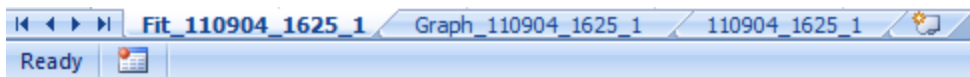

# Fitting functions and parameters

Shirley, Tougaard, Polynomial, Poly Shirley, Poly Tougaard, Victoreen, Arctan

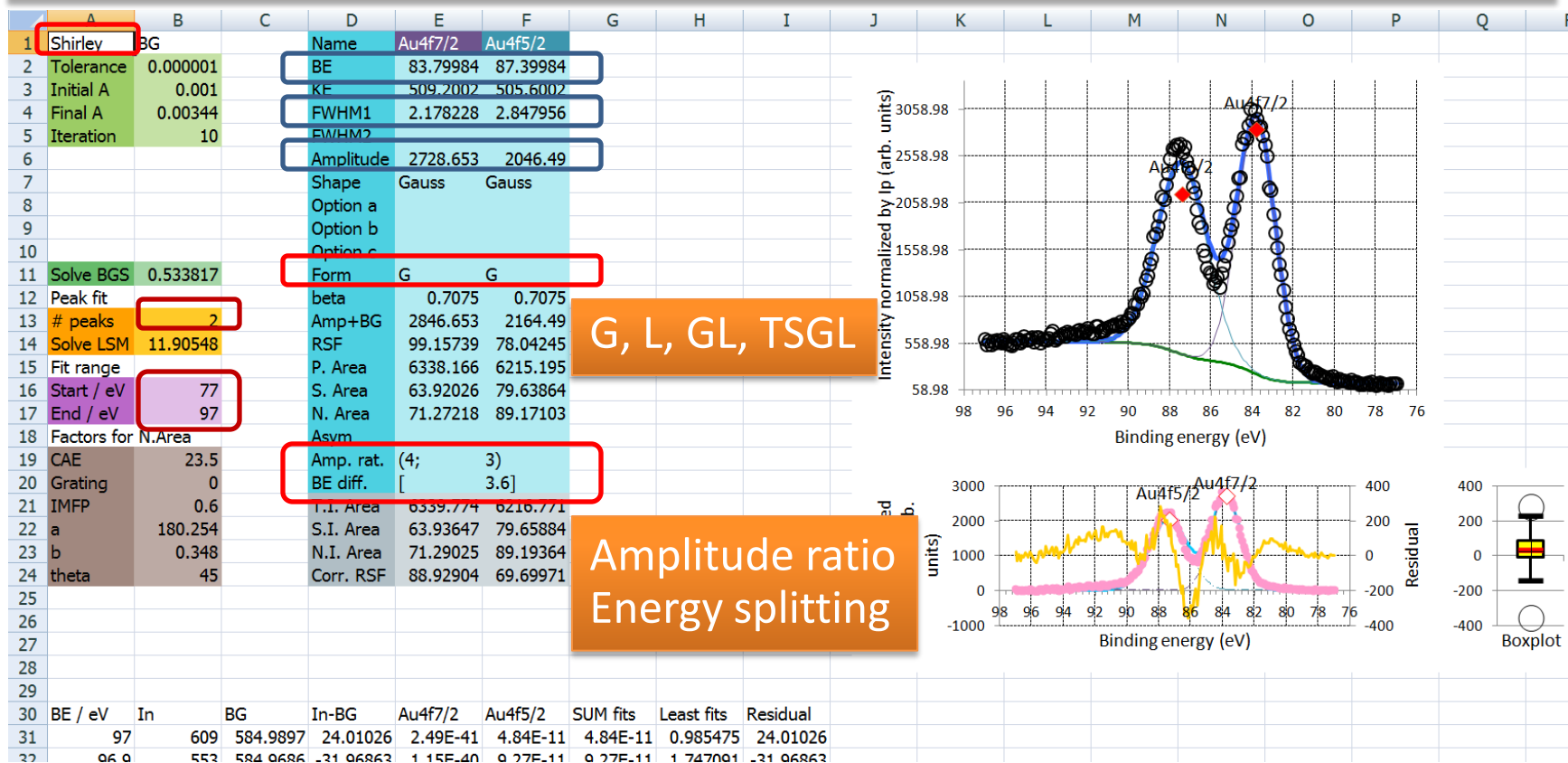



Update results by setup **constraints with bold font style**, various BGs or Forms





# Tailed pseudo-Voigt





	A	B	C	D	E	F	G	H	I	J	K	L	M	N	O	P	Q	R	S
1	Shirley	BG		Name	Au4f7/2	Au4f5/2													
2	Tolerance	0.000001		BE	83.89177	87.49177													
3	Initial A	0.001		KE	509.1082	505.5082													
4	Final A	0.0036		FWHM1	2.192026	2.036631													
5	Iteration	3		FWHM2															
6				Amplitude	88.24273	66.18205													
7				Shape	0.207135	0.001													
8				Option a	0.0927	0.442394													
9				Option b	20.46306	11.55928													
10				Option c															
11	Solve BGS	0.01586		Form	TSGL	TSGL													
12	Peak fit			beta	0.7075	0.7075													
13	# peaks	2		Amp+BG	91.91117	69.85049													
14	Solve LSM	0.045204		RSF	99.15739	78.04245													
15	Fit range			P. Area	226.4484	143.8036													
16	Start / eV	77		S. Area	2.283727	1.842633													
17	End / eV	97		N. Area	2.546665	2.063409													
18	Factors for N.Area			Asym															
19	CAE	23.5		Amp. rat.	(4;	3)													
20	Grating	0		BE diff.	[	3.6]													
21	IMFP	0.6		T.I. Area	240.0941	171.4606													
22	a	180.254		S.I. Area	2.421343	2.197018													
23	b	0.348		N.I. Area	2.700126	2.460255													
24	theta	45		Corr. RSF	88.91957	69.69223													
25																			
26																			
27																			
28																			
29																			
30	BE / eV	In	BG	In-BG	Au4f7/2	Au4f5/2	SUM fits	Least fits	Residual										
31	97	19.69872	18.94799	0.750731	0.271886	0.111915	0.383802	0.006965	0.36693										
32	96.9	17.95129	18.94713	-0.99584	0.278361	0.118642	0.397003	0.100289	-1.39284										
33	96.8	19.18263	18.94387	0.238762	0.285022	0.125775	0.410797	0.001529	-0.17204										
34	96.7	18.03784	18.94352	-0.90568	0.291876	0.133338	0.425214	0.091451	-1.3309										
35	96.6	19.03867	18.94303	0.095641	0.298928	0.141357	0.440284	0.006128	-0.34464										
36	96.5	18.80417	18.94164	-0.13747	0.306184	0.149859	0.456043	0.01816	-0.59352										

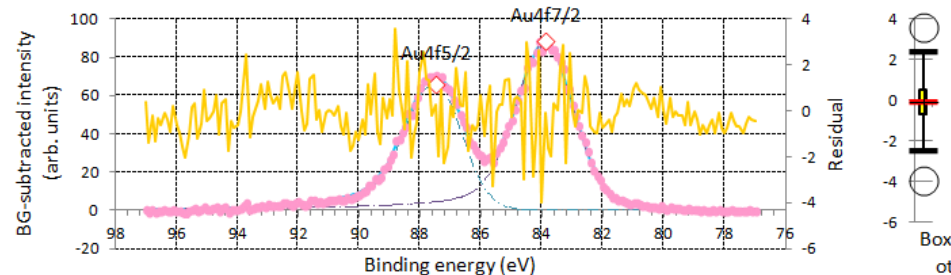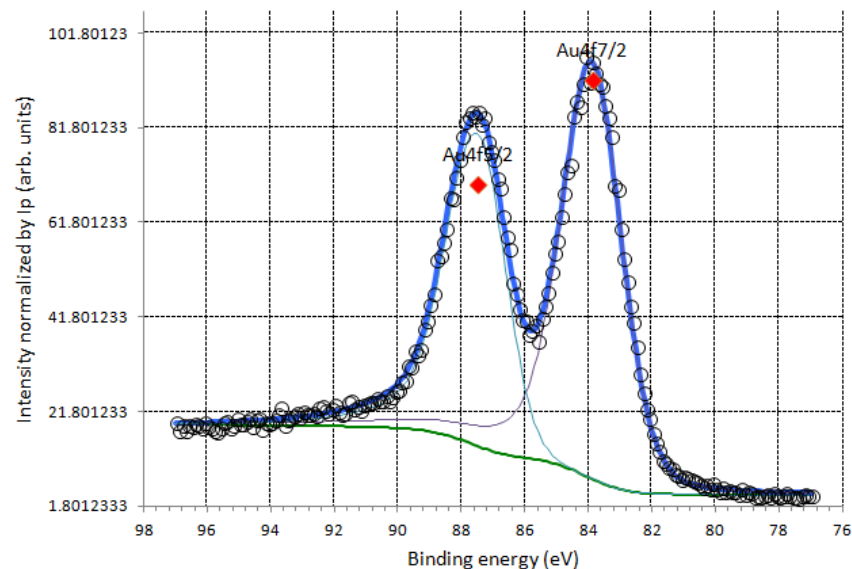


Shirley's BG, Tailed Voigt, Doublet 4f (4:3),  $\Delta$ SO Au4f: 3.6 eV are used for above fitting. User-defined peak and BG functions can be easily implemented in your Visual Basic code.

# Polynomial (Shirley or Tougaard)

## Active-BGs

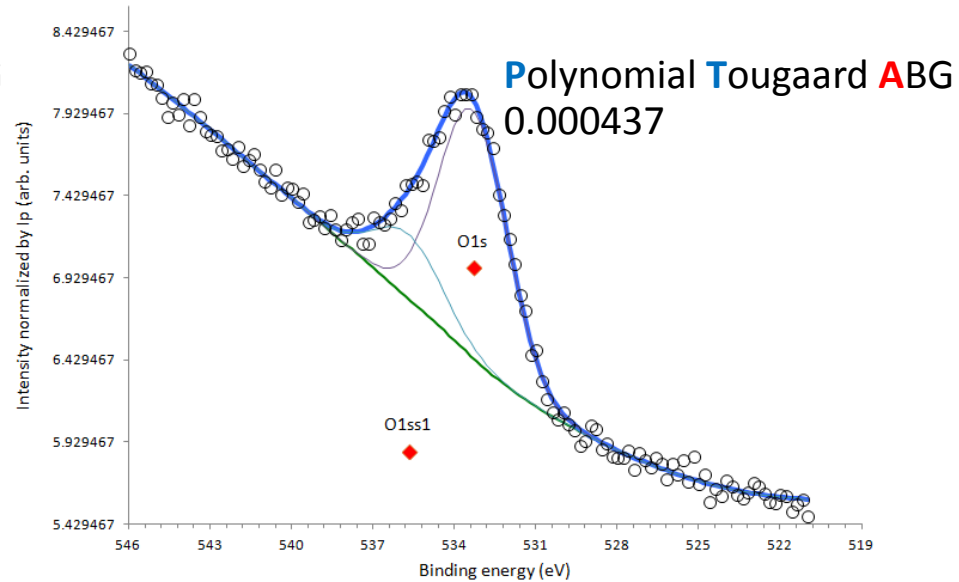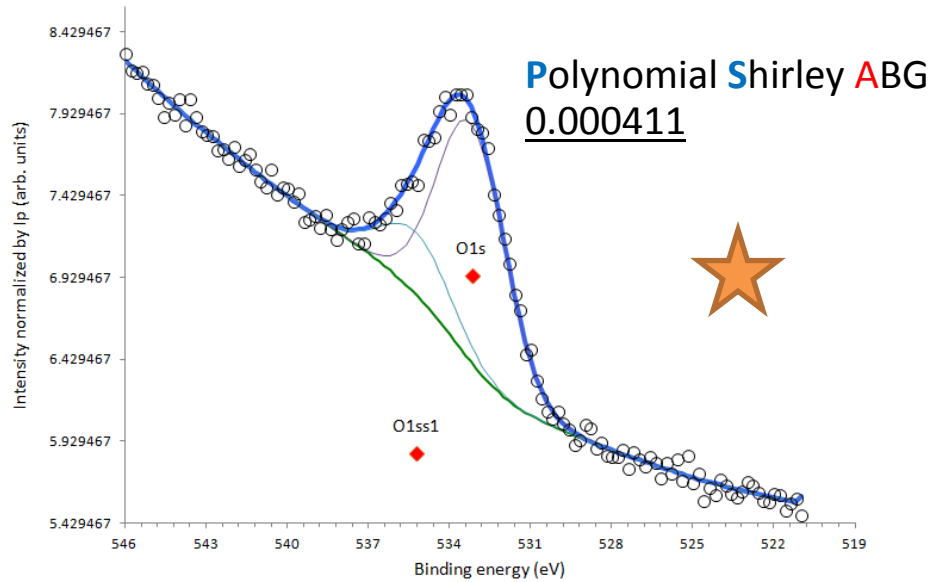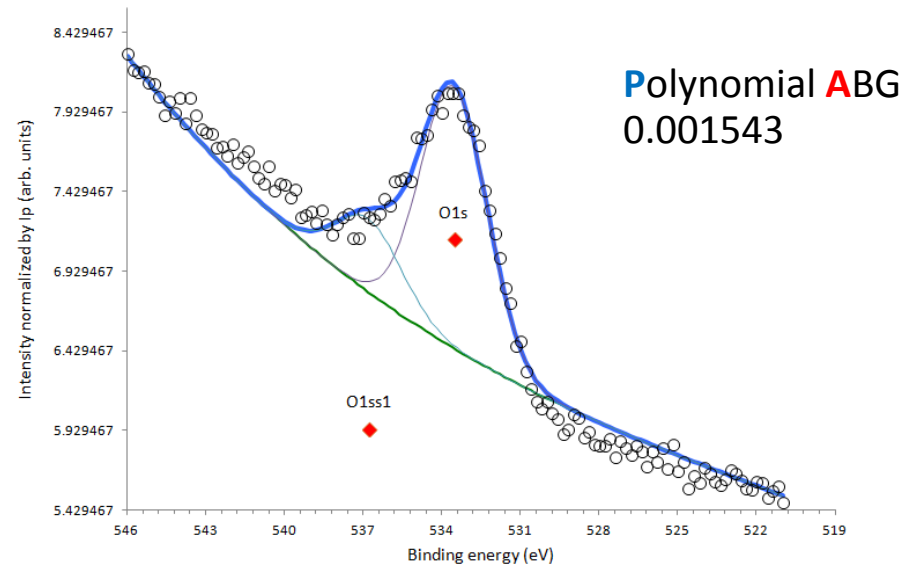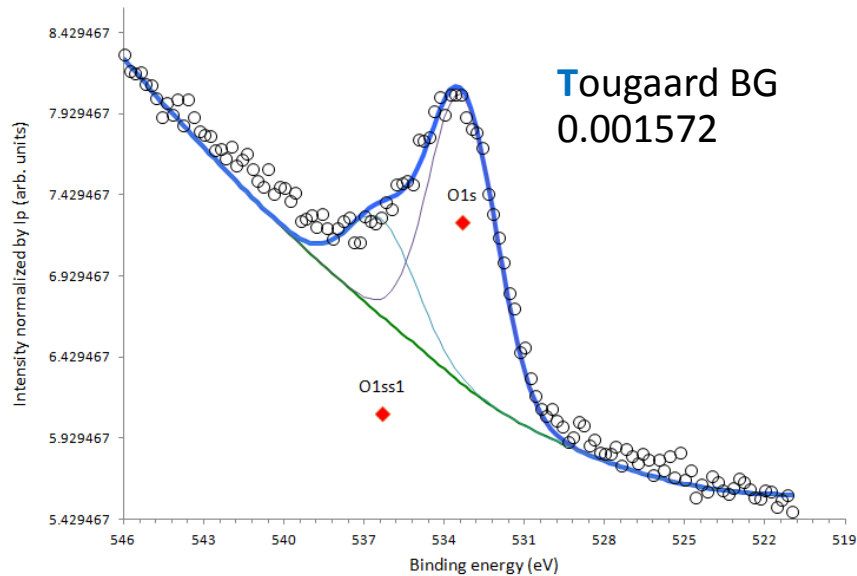

# Poly Shirley ABG

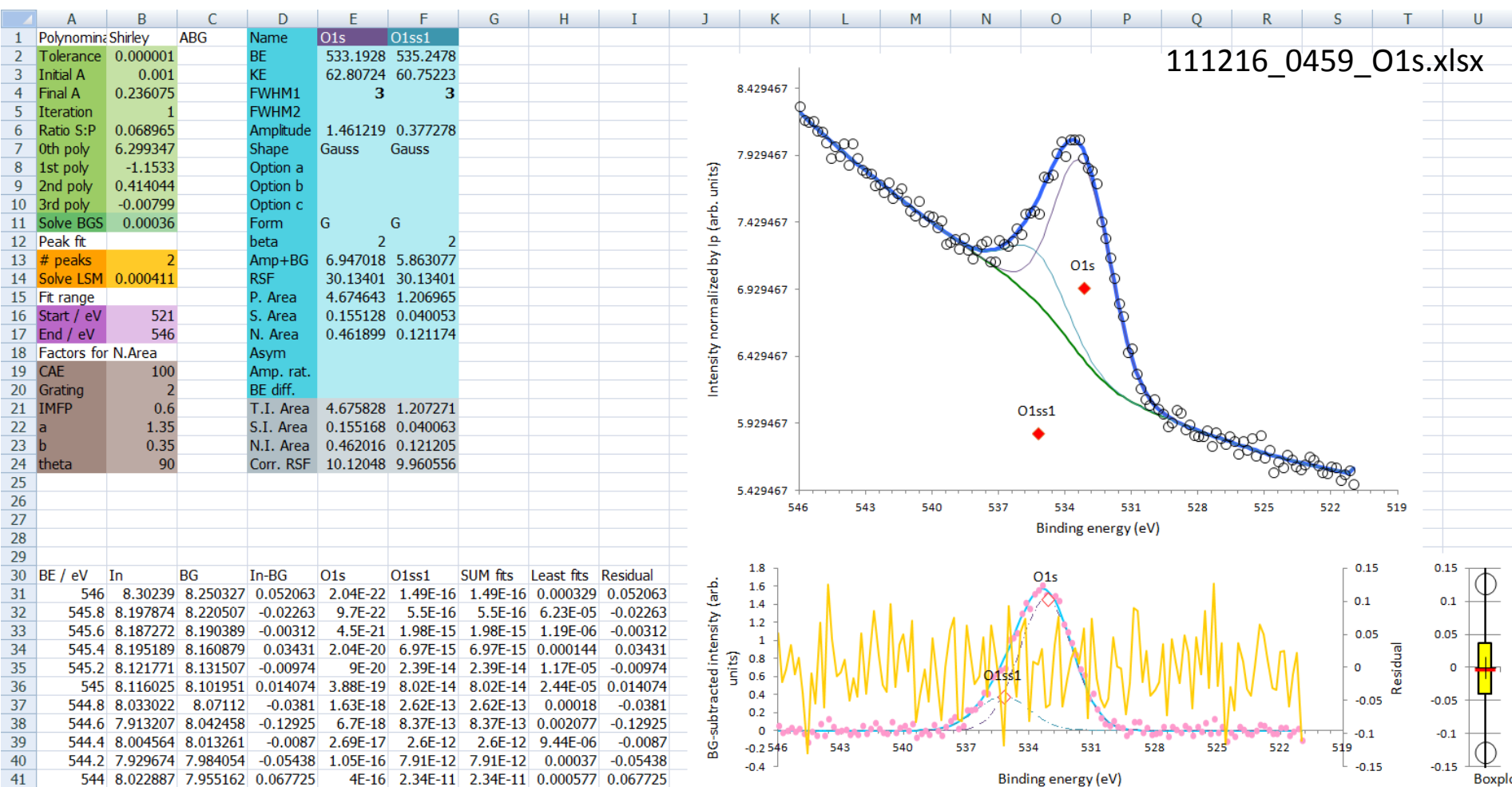

# Lists of fitting functions

## Peak shapes (total 4 forms)

Code (form)	Shape	Option a	Option b
G (0)	Gauss		
L (1)	Lorentz		
<u>TSGL</u>	Exponential blend SGL	Tail scale	Tail length at half max
GL (0 < shape < 1)	G + L with the same FWHM		

## Backgrounds (total 8 BGs)

BG (A1)	BG (B1)	BG (C1)	Technique	Optimization
Shirley	BG		XPS	Static
Tougaard	BG		XPS	Static
Polynomial	BG		XPS/XAS	Static
Polynomial	ABG		XPS/XAS	Active
Polynomial	Shirley	ABG	XPS	Active
Polynomial	Tougaard	ABG	XPS	Active
Victoreen	BG		XAS	Static
ArcTangent	BG		XAS	Active

# Multiple data-fit

	A	B	C	D	E	F
1	Shirley	BG		debug	Au4f7/2	Au4f5/2
2	Tolerance	0.000001		BE	83.89131	87.49132
3	Initial A	0.001		KE	509.1087	505.5087
4	Final A	0.00344		FWHM1	2.196802	2.052602
5	Iteration	6		FWHM2		

“debug”/“debuga” in D1 cell


- debug: apply the original fit parameters
- debuga: apply those used just before




Once fitting process above done, “ana” in D1 cell to compare those fitted xlsx files.




	A	B	C	D	E	F
1	Shirley	BG		ana	Au4f7/2	Au4f5/2
2	Tolerance	0.000001		BE	83.89131	87.49132
3	Initial A	0.001		KE	509.1087	505.5087
4	Final A	0.00344		FWHM1	2.196802	2.052602
5	Iteration	6		FWHM2		

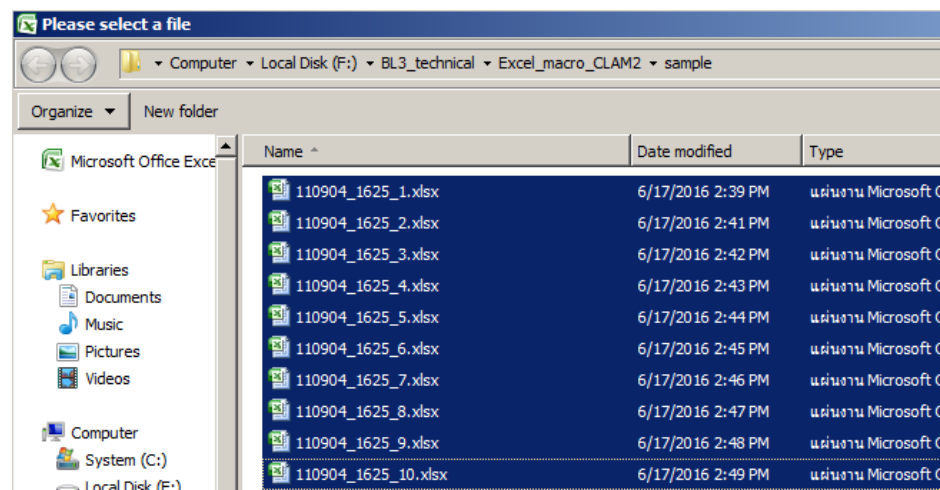

# Summary of fitting results

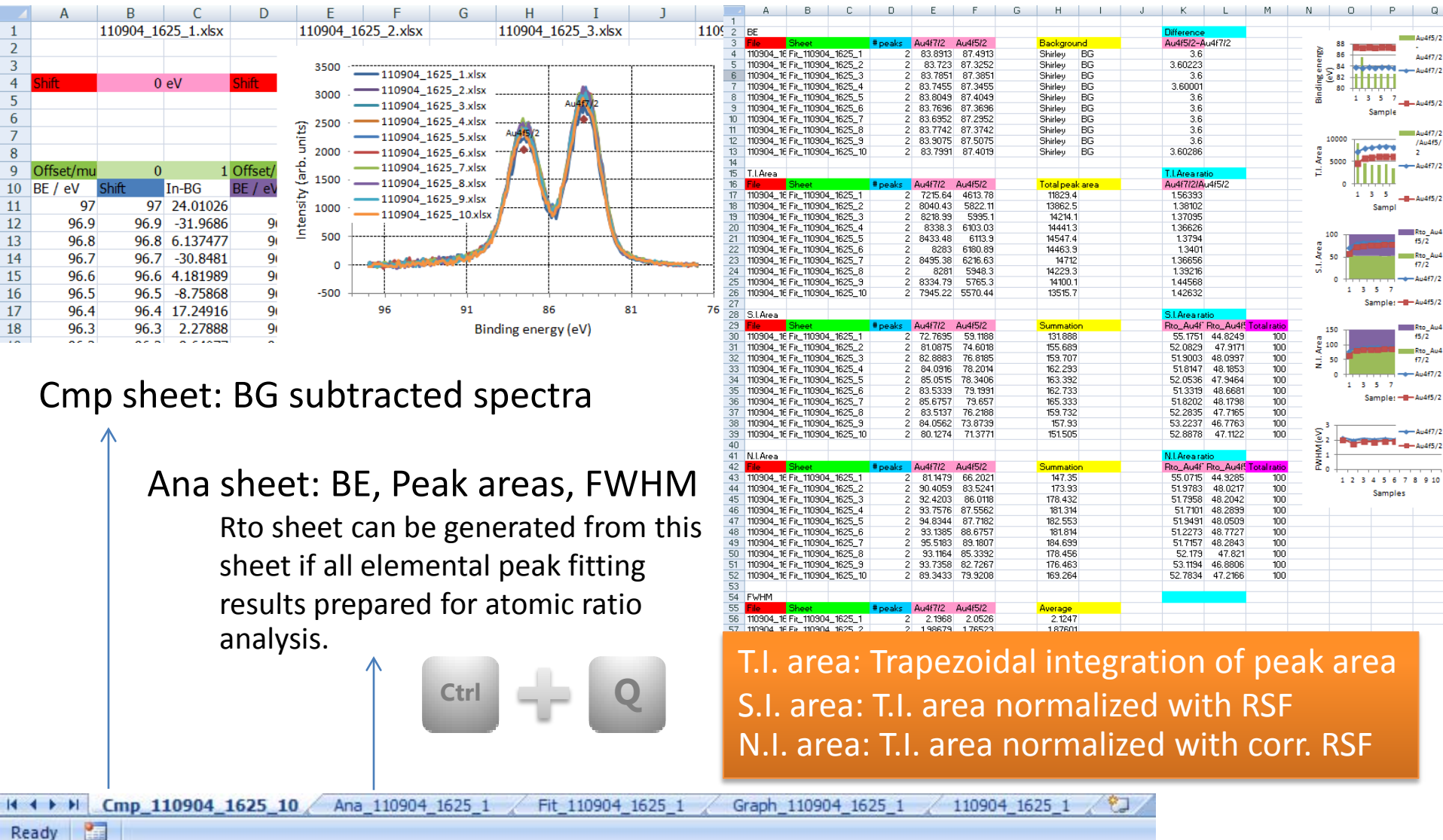

# Enjoy!

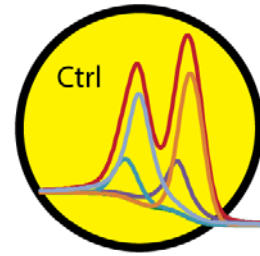

Depth-profile XPS: under review in RSC Adv.

