

Data analysis: Ctrl+Q

- *efficient and effective analysis ever experienced* -

Hideki NAKAJIMA

heitler@gmail.com


18 Aug. 2016


Preface


- **Ctrl+Q** is the first comprehensive data-analysis code for **SR-XPS** and **XAS** working on Windows Office **Excel 2007 or later**.
- Ctrl+Q will analyze the data by just pressing the shortcut key: **Ctrl + q**.
- It works great for the number of publications.
 - Carbon, JPAP, ASS, ACS AMI, SAB etc. at BL3.2Ua in SPL
- **GitHub: xps-excel-macro** for updated info.

 [heitler / xps-excel-macro](#)

 Watch  Star  Fork

 Code

 Issues 0

 Pull requests 0

 Wiki

 Pulse

 Graphs

 Settings


The VBA codes to be pasted in the personal macro workbook (Personal.xlsb) in Windows Office Excel 2007 - 2016 Visual Basic Editor (VBE) work with Solver for spectral analysis and curve fitting of the synchrotron-based soft x-ray photoemission (XPS) and absorption (XAS) spectra.

— Edit

Data preparation

	A	B
1	KE/eV	Scan#1
2	496	19.69872
3	496.1	17.95129
4	496.2	19.18263
5	496.3	18.03784
6	496.4	19.03867
7	496.5	18.80417
8	496.6	18.22572


- Two-column data in the Tab-delimited format texts or spreadsheets are analyzed in the code.
- Sheet is named after the workbook name.
- Syntax in A1 cell has to be chosen as follows.

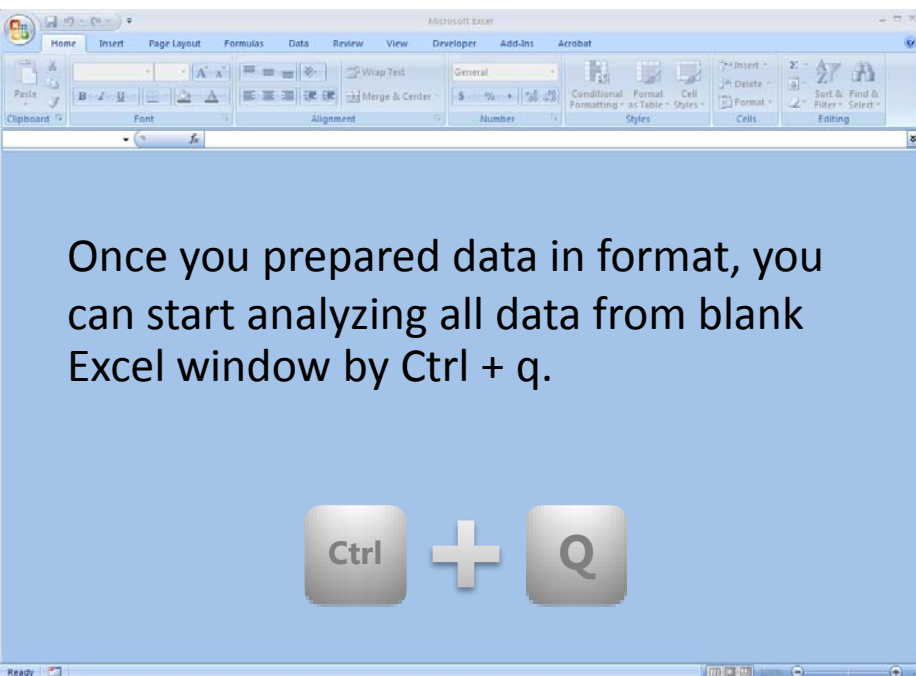

Syntax in A1 cell	Data in column A	Data in column B
KE/eV	KE	XPS
BE/eV	BE	XPS
PE/eV	PE	XAS
GE/eV	PE	G scan
AE/eV	EE	Auger
QE/eV	mass	Q-mass
ME/eV	Position	Any

Installation of the code in Excel

- Generate **Personal Workbook Macro**
- Install **Solver** Add-in and setup its DLL
- Setup the **Shortcut key** in Macro Option

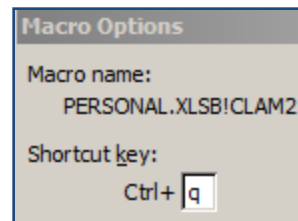

- Compile your own **database** for BE and sensitivities on UD.xlsx


XPS sheet

	A	B	C	D
1	Element	Orbit	BE(eV)	ASF
2	C	1s	284.6	1
3	O	1s	532	2.93

AES sheet

	A	B	C	D
1	Element	Auger	KE(eV)	RSF
2	C	KLL	266	0.6
3	O	KLL	506	0.96


C1s sample csv

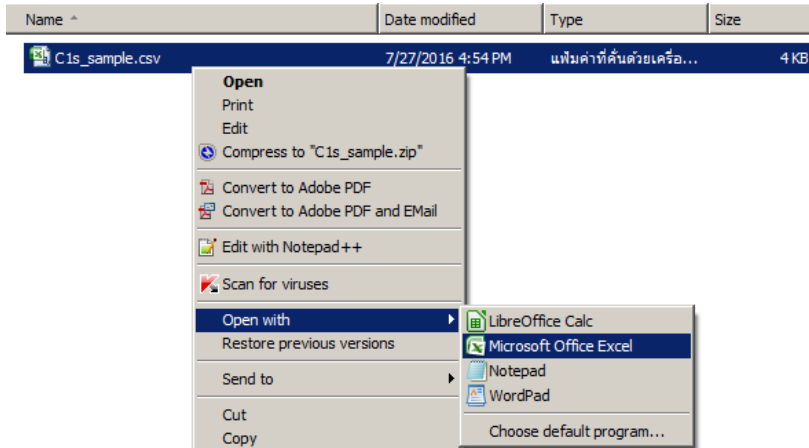



	A	B
1	KE/eV	PE: 600 eV
2	290	1709
3	290.1	1798
4	290.2	1845
5	290.3	1722
6	290.4	1729
7	290.5	1699
8	290.6	1727
9	290.7	1693
10	290.8	1647
11	290.9	1684
12	291	1650



600 eV of photon energy

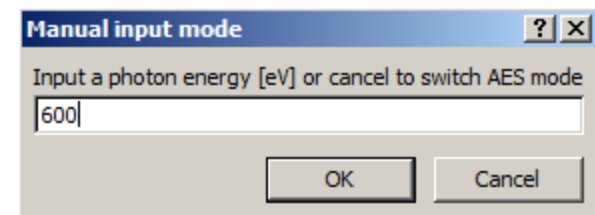

XPS: Carbon 1s spectrum

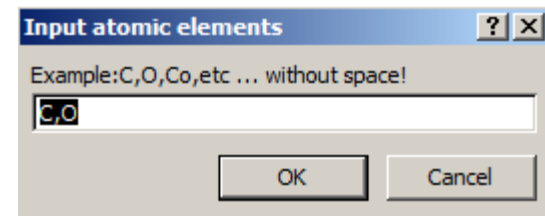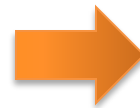

Sample file with multiple scans

File Explorer context menu for **Au4f_sample.csv**:

- Open
- Print
- Edit
- Compress to "Au4f_sample.zip"
- Convert to Adobe PDF
- Convert to Adobe PDF and Email
- Edit with Notepad ++
- Scan for viruses
- Open with
 - LibreOffice Calc
 - Microsoft Office Excel**
 - Notepad
 - WordPad
- Restore previous versions
- Send to
 - Choose default program...
- Cut
- Copy
- Create shortcut
- Delete
- Rename
- Properties

Microsoft Excel screenshot showing a table with columns A through H and rows 1 through 7. The table contains energy values (KE/eV) and some cells are marked with 'X' and 'Y'.

File Explorer screenshot showing a list of files named 110904_1625_1.txt through 110904_1625_11.txt.

Microsoft Excel window showing a warning message: "Data were exported in the text files."

- Download sample.txt including 100 spectra.
- Open with Excel as shown above figure.
- Export 100 text files by shortcut keys Ctrl+Q.

Sample file with single scan

Name ^	Date modified	Type	Size
110904_1625_1.txt	7/27/2016 5:12 PM	TXT File	3 KB
110904_1625_2.txt		TXT File	3 KB
110904_1625_3.txt		TXT File	3 KB
110904_1625_4.txt		TXT File	3 KB
110904_1625_5.txt		TXT File	3 KB
110904_1625_6.txt		TXT File	3 KB
110904_1625_7.txt		TXT File	3 KB
110904_1625_8.txt		TXT File	3 KB
110904_1625_9.txt		TXT File	3 KB
110904_1625_10.txt		TXT File	3 KB
110904_1625_11.txt		TXT File	3 KB
110904_1625_12.txt		TXT File	3 KB
110904_1625_13.txt		TXT File	3 KB
110904_1625_14.txt		TXT File	3 KB
110904_1625_15.txt		TXT File	3 KB
110904_1625_16.txt		TXT File	3 KB
110904_1625_17.txt		TXT File	3 KB
110904_1625_18.txt		TXT File	3 KB

Open
Edit
Compress to "110904_1625_1.zip"
Convert to Adobe PDF
Combine supported files in Acrobat...
Edit with Notepad++
Scan for viruses

Open with ▶
Restore previous versions

Send to ▶


Cut
Copy

Create shortcut
Delete
Rename

Properties

Adobe Illustrator CS6
LibreOffice Calc
LibreOffice Writer
Mery
Microsoft Office Excel
Notepad
Notepad++ : a free (GNU) source co
WordPad

Choose default program...

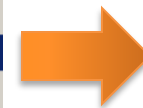

	A	B
1	KE/eV	110904_16
2	496	609
3	496.1	553
4	496.2	591
5	496.3	554
6	496.4	589
7	496.5	576
8	496.6	602
9	496.7	587
10	496.8	576
11	496.9	608
12	497	594

Data to plot in Graph sheet

	A	B
1	KE/eV	110904.16
2	496	609
3	496.1	553
4	496.2	591
5	496.3	554
6	496.4	589
7	496.5	576
8	496.6	602
9	496.7	587
10	496.8	576
11	496.9	608
12	497	584

KE/eV represents a kinetic energy.



600 eV of photon energy →

Manual input mode ? X

Input a photon energy [eV] or cancel to switch AES mode

600

OK Cancel

Sample: Gold metal foil →

Input atomic elements ? X

Example: C,O,Co,etc ... without space!

Au

OK Cancel

Note that Au4f sensitivity is not in database, so you have to input factors in database in a way below.

Au 4f BE and ASF
in XPS sheet on UD.xlsx →

	A	B	C	D
1	Element	Orbit	BE(eV)	ASF
2	C	1s	284.6	1
3	O	1s	532	2.93
4	Au	4f5/2	87.6	7.54
5	Au	4f7/2	84	9.58

[illegible]


Ctrl + Q

Offset/multiplied data in bottom

Reference spectrum

Grating	0 lines/mm	
PE	600 eV	
WF	4 eV	
Char	3 eV	
Start KE	496 eV	
End KE	516 eV	
Step KE	0.1 eV	
# scan	1 times	
Offset/mu	0	1
Ke	Be	In
496	97	609
496.1	96.9	553
496.2	96.8	591
496.3	96.7	554
496.4	96.6	589
496.5	96.5	576
496.6	96.4	602
496.7	96.3	587
496.8	96.2	576

Standard peak BE should be well-known.

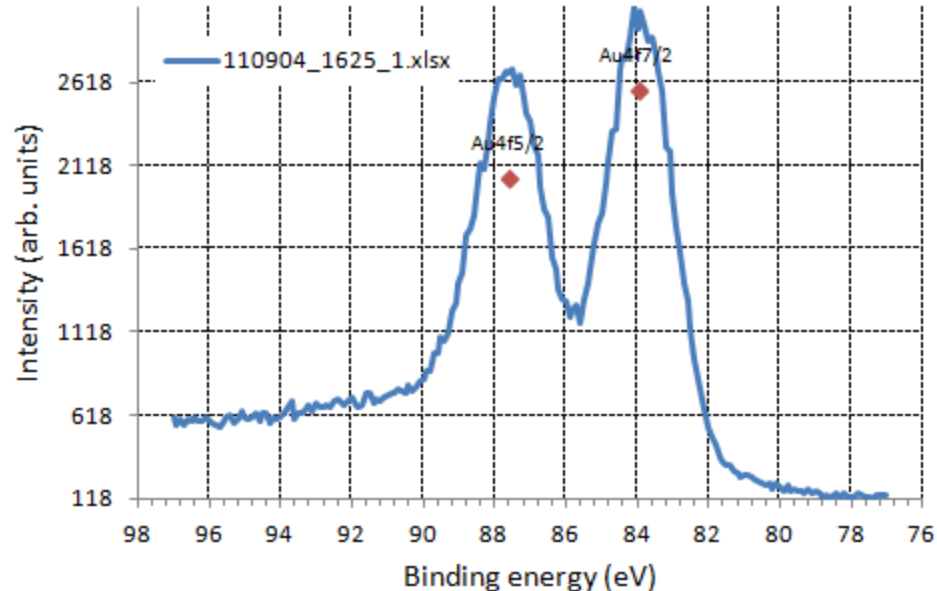


Next, apply these factors to another by
“debug” in A1 cell then



	A	B	C
1	debug	0 lines/mm	
2	PE	600 eV	
3	WF	4 eV	
4	Char	3 eV	
5	Start KE	496 eV	
6	End KE	516 eV	
7	Step KE	0.1 eV	
8	# scan	1 times	
9	Offset/mu	0	1
10	Ke	Be	In

Multiple-file selection

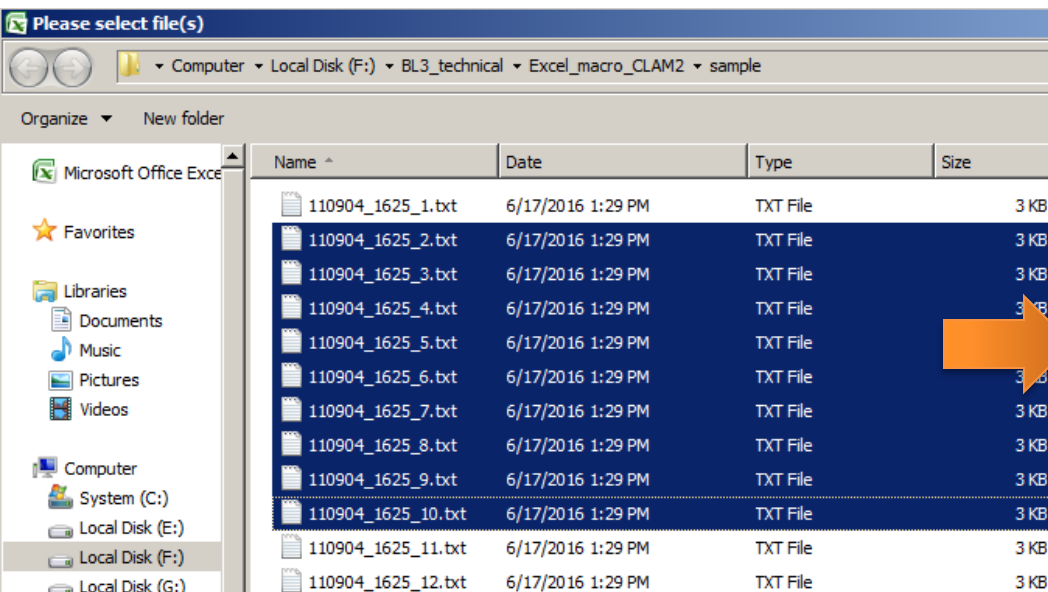

Name ^	Date	Type	Size
110904_1625_1.txt	6/17/2016 1:29 PM	TXT File	3 KB
110904_1625_1.xlsx	6/17/2016 1:55 PM	แผ่นงาน Microsoft O...	53 KB
110904_1625_2.txt	6/17/2016 1:29 PM	TXT File	3 KB
110904_1625_2.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	52 KB
110904_1625_3.txt	6/17/2016 1:29 PM	TXT File	3 KB
110904_1625_3.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	12 KB
110904_1625_4.txt	6/17/2016 1:29 PM	TXT File	3 KB
110904_1625_4.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	53 KB
110904_1625_5.txt	6/17/2016 1:29 PM	TXT File	3 KB
110904_1625_5.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	53 KB
110904_1625_6.txt	6/17/2016 1:29 PM	TXT File	3 KB
110904_1625_6.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	52 KB
110904_1625_7.txt	6/17/2016 1:29 PM	TXT File	3 KB
110904_1625_7.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	52 KB
110904_1625_8.txt	6/17/2016 1:29 PM	TXT File	3 KB

Next, compare the spectra by “comp” in D1

	A	B	C	D
1	Grating	0 lines/mm	comp	
2	PE	600 eV		
3	WF	4 eV		
4	Char	3 eV		
5	Start KE	496 eV		
6	End KE	516 eV		
7	Step KE	0.1 eV		
8	# scan	1 times		
9	Offset/mu	0	1	
10	Ke	Be	In	
11	496	97	609	
12	496.1	96.9	553	
13	496.2	96.8	591	
14	496.3	96.7	554	
15	496.4	96.6	589	
16	496.5	96.5	576	
17	496.6	96.4	602	
18	496.7	96.3	587	




File Explorer window showing a directory of files. The 'File name' column shows files named 110904_1625_1.xlsx through 110904_1625_10.xlsx. The 'Date modified' column shows dates from 6/17/2016. The 'Type' column shows 'แผ่นงาน Microsoft O...'. The 'Size' column shows various sizes from 53 KB to 181 KB. An orange arrow points from the file list to the right.

Name ^	Date modified	Type	Size
110904_1625_1.xlsx	6/17/2016 2:06 PM	แผ่นงาน Microsoft O...	177 KB
110904_1625_2.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	53 KB
110904_1625_3.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	53 KB
110904_1625_4.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	53 KB
110904_1625_5.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	53 KB
110904_1625_6.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	52 KB
110904_1625_7.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	52 KB
110904_1625_8.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	52 KB
110904_1625_9.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	52 KB
110904_1625_10.xlsx	6/17/2016 2:02 PM	แผ่นงาน Microsoft O...	53 KB
sample.xlsx	6/17/2016 1:55 PM	แผ่นงาน Microsoft O...	181 KB

Compared spectra

Grating	0 lines/mm		110904_1625_2.xlsx			110904_1625_3.xlsx			1109
PE	600 eV		PE	600 eV		PE	600 eV		PE
WF	4 eV		WF	4 eV		WF	4 eV		WF
Char	3 eV		Char	3 eV		Char	3 eV		Char
Start KE	496 eV		Start KE	496 eV		Start KE	496 eV		Start KE
End KE	516 eV		End KE	516 eV		End KE	516 eV		End KE
Step KE	0.1 eV		Step KE	0.1 eV		Step KE	0.1 eV		Step KE
# scan	1 times		# scan	1 times		# scan	1 times		# scan
Offset/mu	0		Offset/mu	0		Offset/mu	0		Offset/mu
Ke	Be	In	Ke	Be	In	Ke	Be	In	Ke

496	97	609
496.1	96.9	553
496.2	96.8	591
496.3	96.7	554
496.4	96.6	589
496.5	96.5	576
496.6	96.4	602
496.7	96.3	587
496.8	96.2	576
496.9	96.1	608
497	96	584
497.1	95.9	566
497.2	95.8	555
497.3	95.7	543
497.4	95.6	576
497.5	95.5	606
497.6	95.4	614
497.7	95.3	575

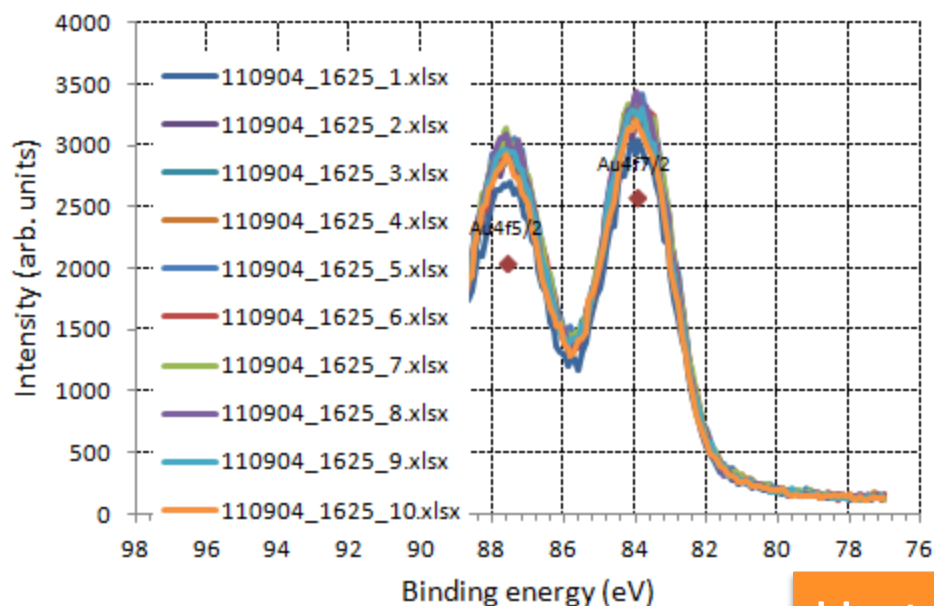

Up to 100 spectra
for BE&KE plots

Spectral comparison functions


- Automatically scaled spectra for comparison
 - “auto” in A1 cell to scale at the both ends
 - “auto[x1:x2,x3:x4]” in A1 cell to scale in the ranges of (x1, x2) for offset (0) and (x3, x4) for multiple (1).
- Annotate legends in each plot
 - “leg” in A1 cell to generate Sample sheet
 - Describe sample or spectrum name in the sheet
 - “leg” in A1 cell to run the code again.
- Normalize spectra with ref. spectrum by “norm” in A1 cell.
 - Second set of data will be selected for ref. and third set to be normalized data.
 - Norm sheet appears from the data normalized.

Every single step needs




	A	B	C		A	B	C	D	E
1	leg	0	lines/mm		1	No.	Name	Sep.	File name
2	PE	600 eV			2	1	Au1	:	110904_1625_1.xlsx
3	WF	4 eV			3	2	Au2	:	110904_1625_2.xlsx
4	Char	3 eV			4	3	Au3	:	110904_1625_3.xlsx
5	Start KE	496			5	4	Au4	:	110904_1625_4.xlsx
6	End KE	516			6	5	Au5	:	110904_1625_5.xlsx
7	Step KE	0.1 eV			7	6	Au6	:	110904_1625_6.xlsx
8	# scan	1 times			8	7	Au7	:	110904_1625_7.xlsx
9	Offset/mu	0	1		9	8	Au8	:	110904_1625_8.xlsx
10	Ke	Be	In		10	9	Au9	:	110904_1625_9.xlsx
11					11	10	Au10	:	110904_1625_10.xlsx


Sample sheet

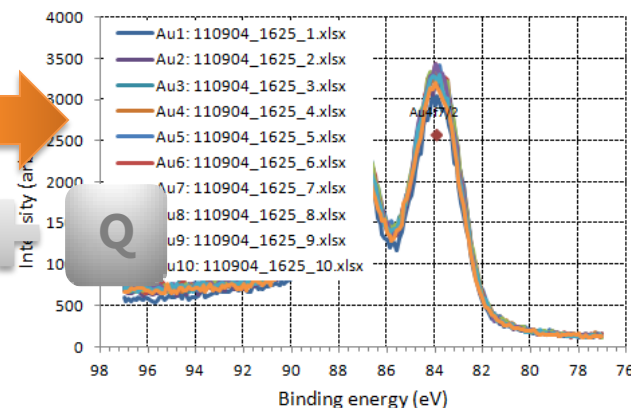

Fit sheet

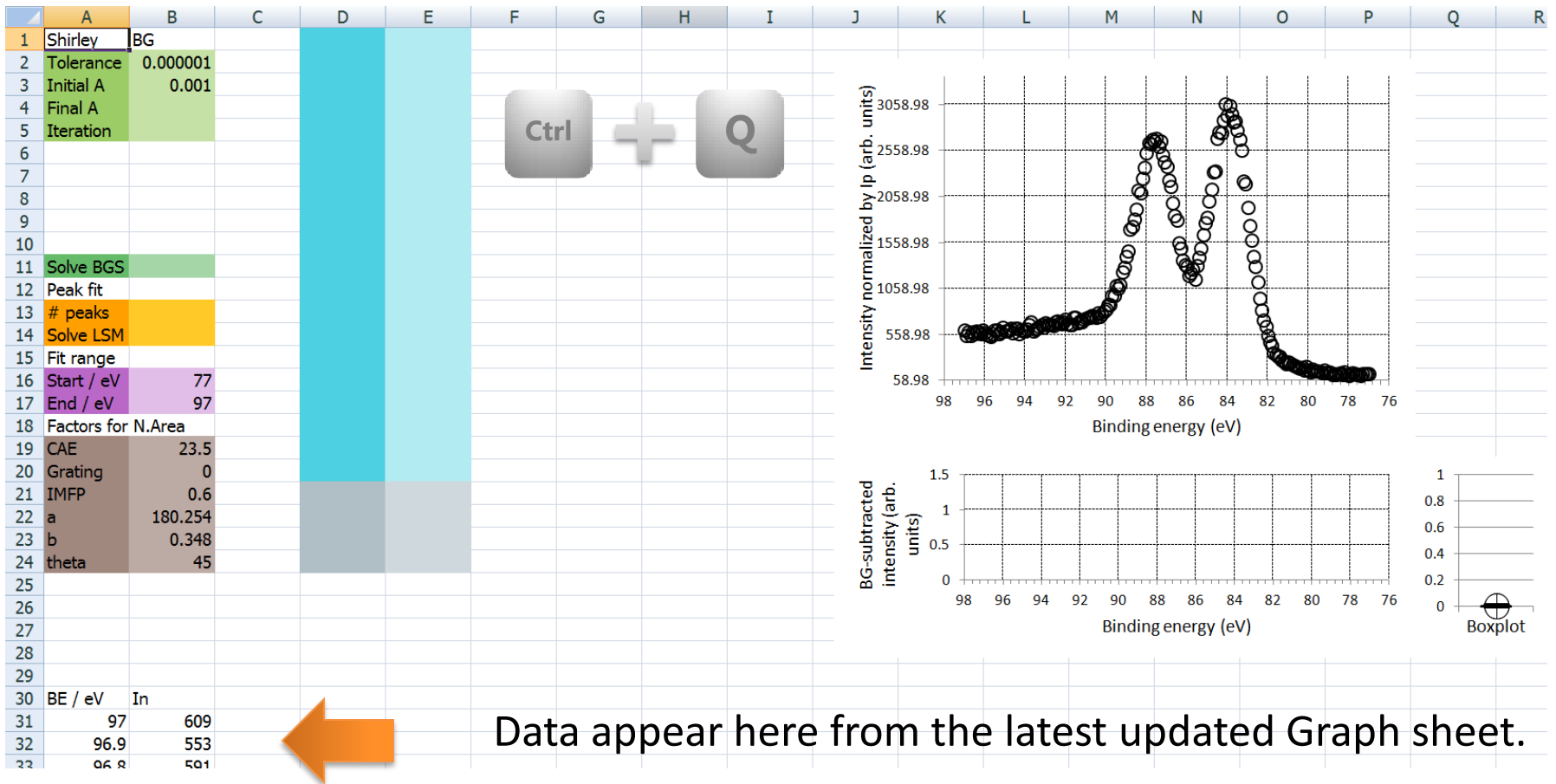


Switch the sheet

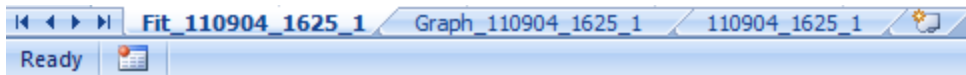

Fitting parameters

S, T, P, PS, PT, V, A, etc.


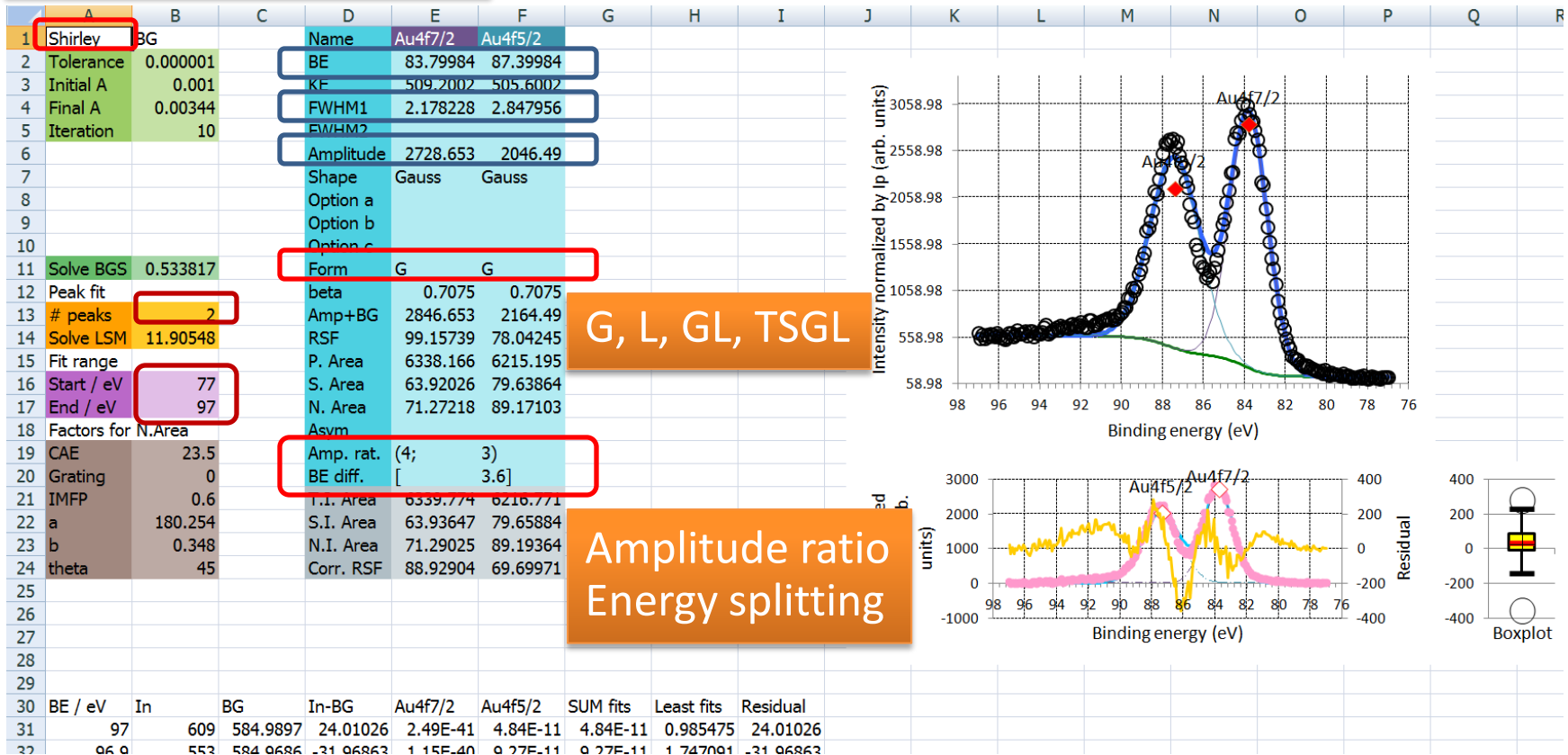


Update results by setup **constraints with bold font style**, various BGs or Forms







Fitting result: tailed pseudo-Voigt

	A	B	C	D	E	F	G	H	I	J	K	L	M	N	O	P	Q	R	S
1	Shirley	BG		Name	Au4f7/2	Au4f5/2													
2	Tolerance	0.000001		BE	83.89177	87.49177													
3	Initial A	0.001		KE	509.1082	505.5082													
4	Final A	0.0036		FWHM1	2.192026	2.036631													
5	Iteration	3		FWHM2															
6				Amplitude	88.24273	66.18205													
7				Shape	0.207135	0.001													
8				Option a	0.0927	0.442394													
9				Option b	20.46306	11.55928													
10				Option c															
11	Solve BGS	0.01586		Form	TSGL	TSGL													
12	Peak fit			beta	0.7075	0.7075													
13	# peaks	2		Amp+BG	91.91117	69.85049													
14	Solve LSM	0.045204		RSF	99.15739	78.04245													
15	Fit range			P. Area	226.4484	143.8036													
16	Start / eV	77		S. Area	2.283727	1.842633													
17	End / eV	97		N. Area	2.546665	2.063409													
18	Factors for N.Area			Asym															
19	CAE	23.5		Amp. rat.	(4;	3)													
20	Grating	0		BE diff.	[3.6]													
21	IMFP	0.6		T.I. Area	240.0941	171.4606													
22	a	180.254		S.I. Area	2.421343	2.197018													
23	b	0.348		N.I. Area	2.700126	2.460255													
24	theta	45		Corr. RSF	88.91957	69.69223													
25																			
26																			
27																			
28																			
29																			
30	BE / eV	In	BG	In-BG	Au4f7/2	Au4f5/2	SUM fits	Least fits	Residual										
31	97	19.69872	18.94799	0.750731	0.271886	0.111915	0.383802	0.006965	0.36693										
32	96.9	17.95129	18.94713	-0.99584	0.278361	0.118642	0.397003	0.100289	-1.39284										
33	96.8	19.18263	18.94387	0.238762	0.285022	0.125775	0.410797	0.001529	-0.17204										
34	96.7	18.03784	18.94352	-0.90568	0.291876	0.133338	0.425214	0.091451	-1.3309										
35	96.6	19.03867	18.94303	0.095641	0.298928	0.141357	0.440284	0.006128	-0.34464										
36	96.5	18.80417	18.94164	-0.13747	0.306184	0.149859	0.456043	0.01816	-0.59352										

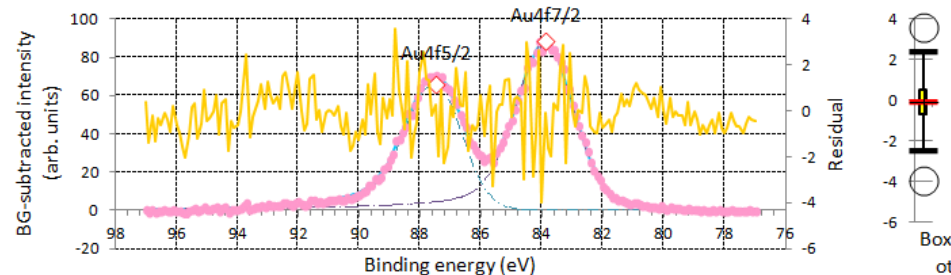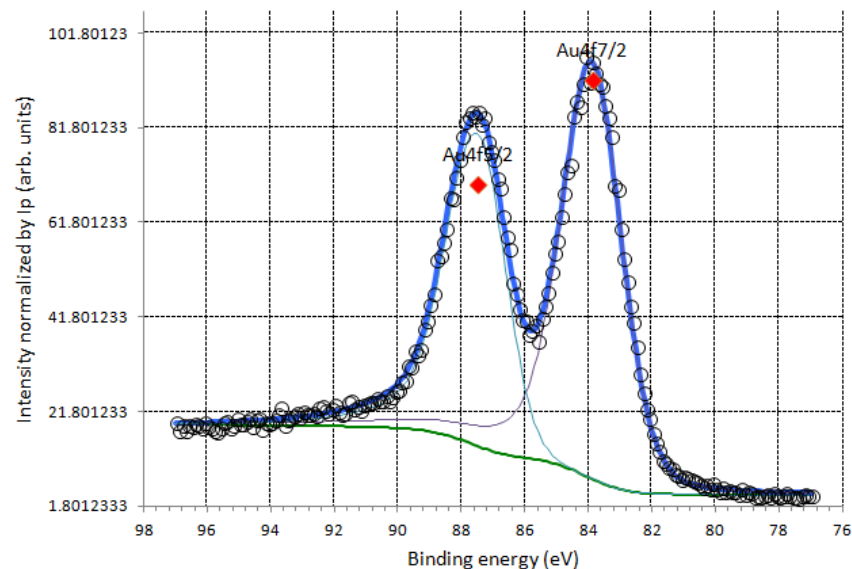


Shirley's BG, Tailed Voigt, Doublet 4f (4:3), ΔSO Au4f: 3.6 eV are used for above fitting. User-defined peak and BG functions can be easily implemented in your Visual Basic code.

Fitting result: Polynomial Shirley & Tougaard Active-BGs

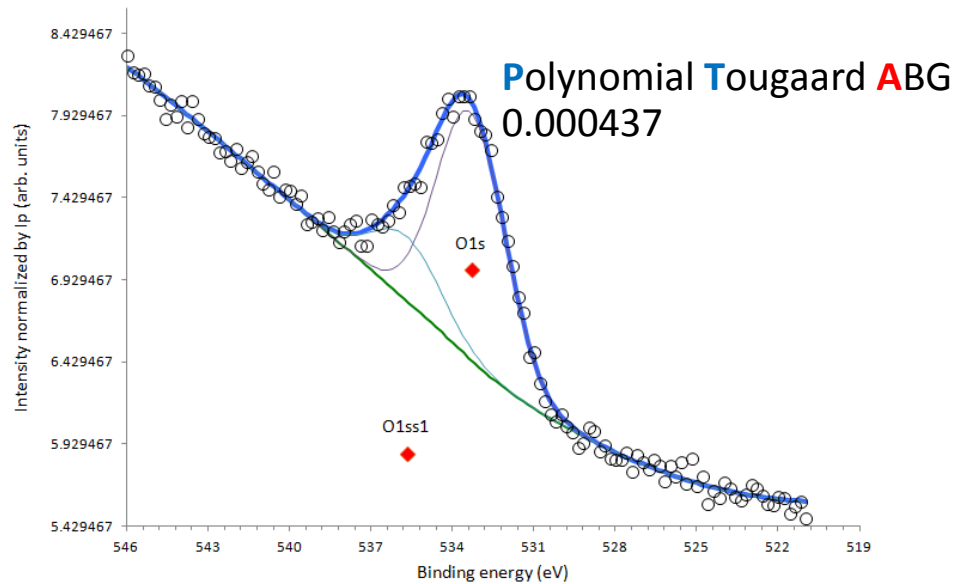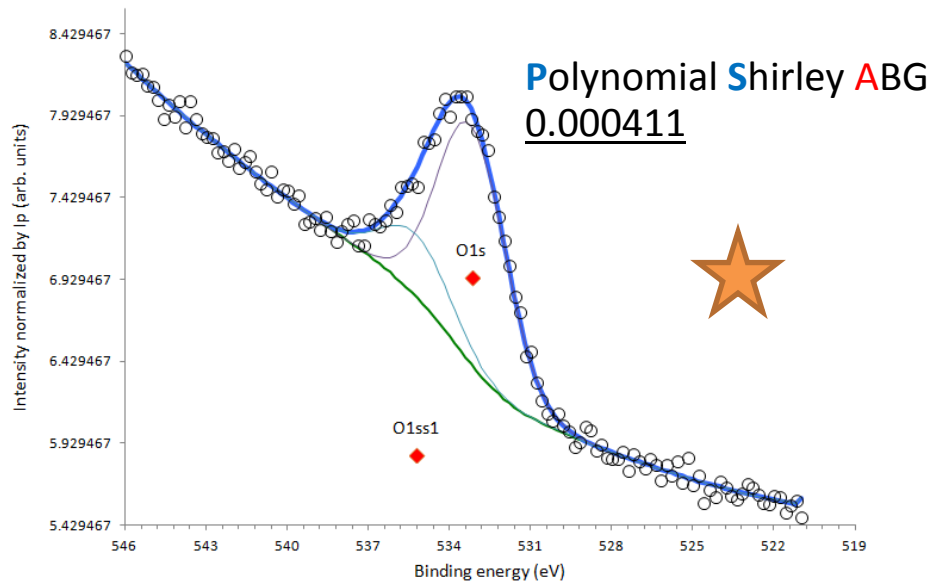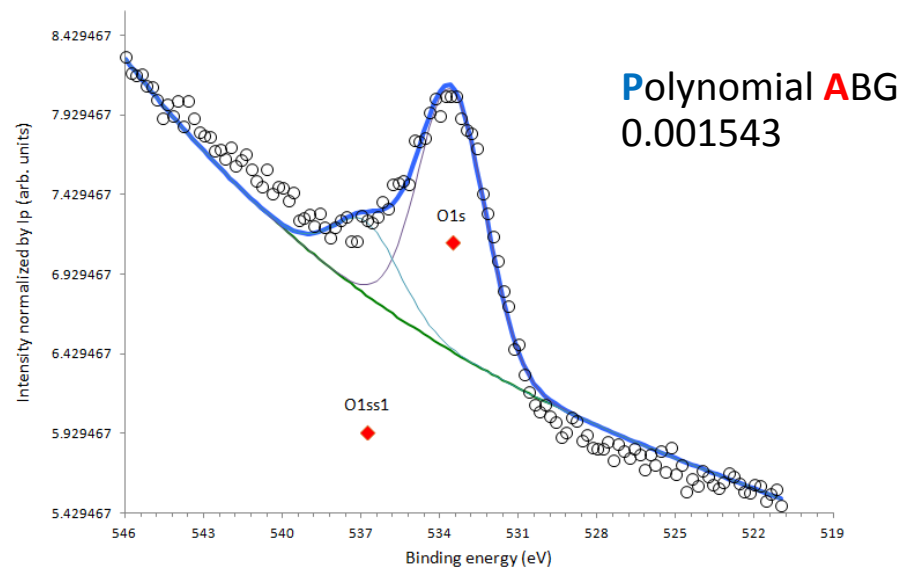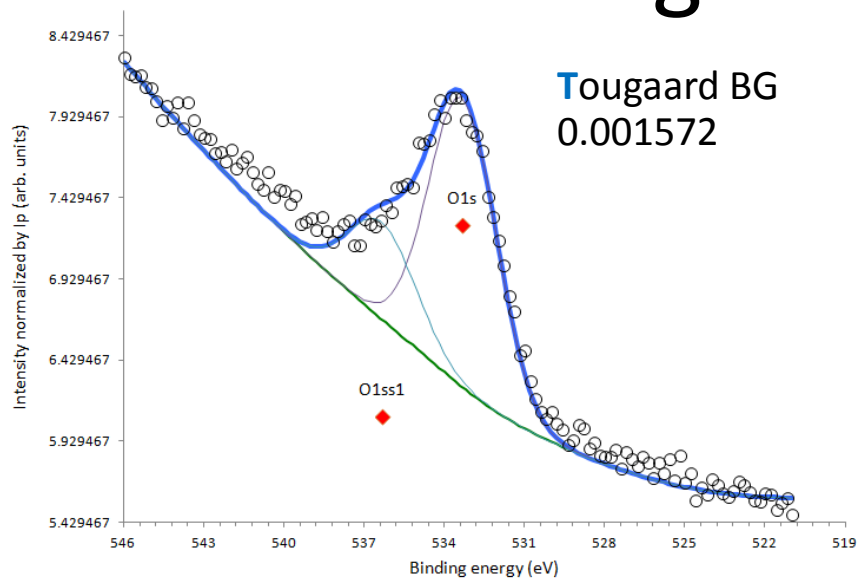

Fitting result: Poly Shirley ABG

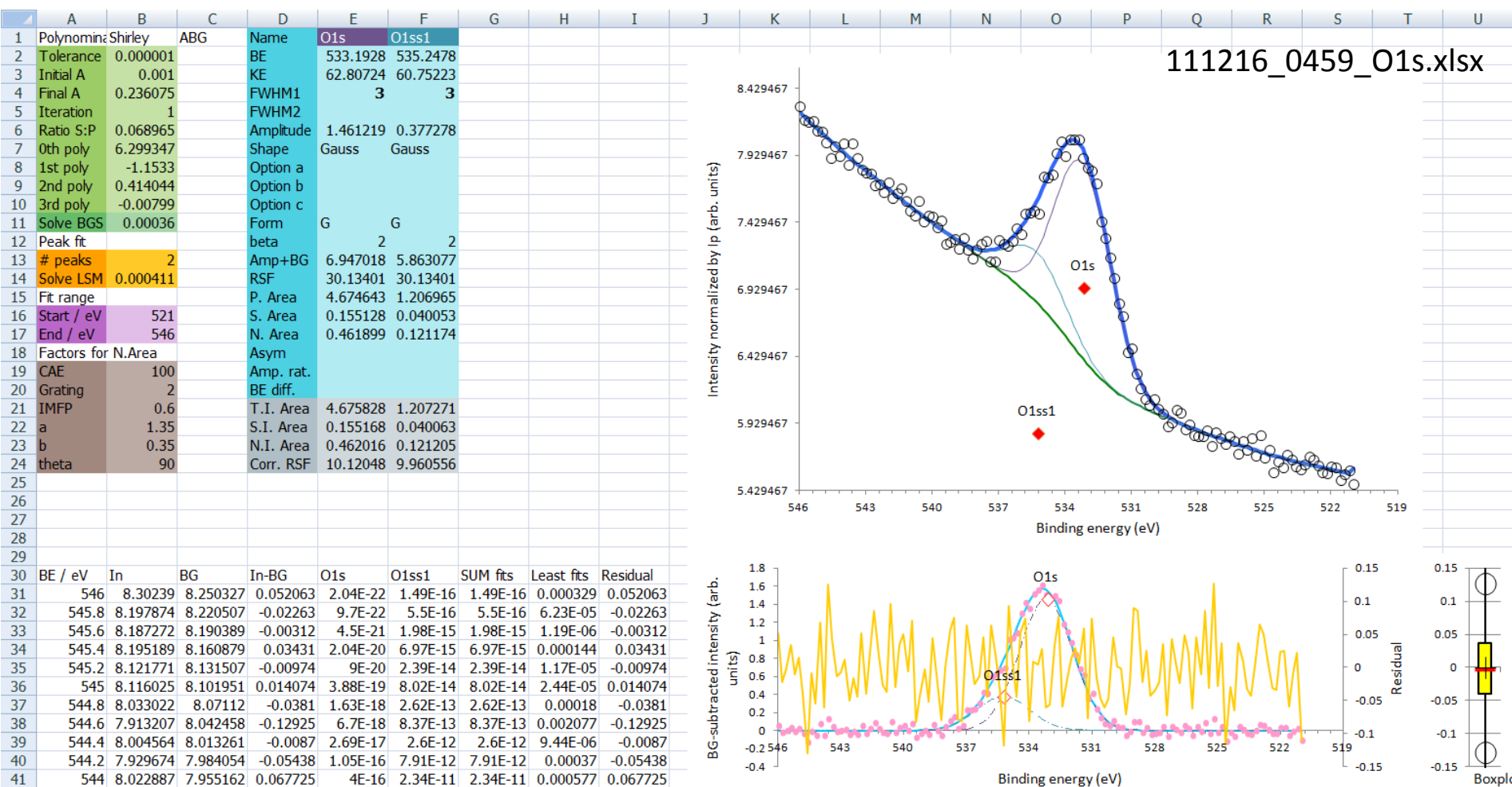

Peak shapes (total 4 forms)

Code (form)	Shape	Option a	Option b
G (0)	Gauss		
L (1)	Lorentz		
<u>TSGL</u>	Exponential blend SGL (MultiPak)	Tail scale	Tail length at half max
GL (0 < shape < 1)	G + L with the same FWHM (MultiPak)		

Backgrounds (total 8 BGs)

BG (A1)	BG (B1)	BG (C1)	Technique	Optimization
Shirley	BG		XPS	Static
Tougaard	BG		XPS	Static
Polynomial	BG		XPS/XAS	Static
Polynomial	ABG		XPS/XAS	Active
Polynomial	Shirley	ABG	XPS	Active
Polynomial	Tougaard	ABG	XPS	Active
Victoreen	BG		XAS	Static
ArcTangent	BG		XAS	Active

Multiple data-fit

	A	B	C	D	E	F
1	Shirley	BG		debug	Au4f7/2	Au4f5/2
2	Tolerance	0.000001		BE	83.89131	87.49132
3	Initial A	0.001		KE	509.1087	505.5087
4	Final A	0.00344		FWHM1	2.196802	2.052602
5	Iteration	6		FWHM2		

“debug”/“debuga” in D1 cell


- debug: apply the original fit parameters
- debuga: apply those used just before




Once fitting process above done, “ana” in D1 cell to compare those fitted xlsx files.




	A	B	C	D	E	F
1	Shirley	BG		ana	Au4f7/2	Au4f5/2
2	Tolerance	0.000001		BE	83.89131	87.49132
3	Initial A	0.001		KE	509.1087	505.5087
4	Final A	0.00344		FWHM1	2.196802	2.052602
5	Iteration	6		FWHM2		

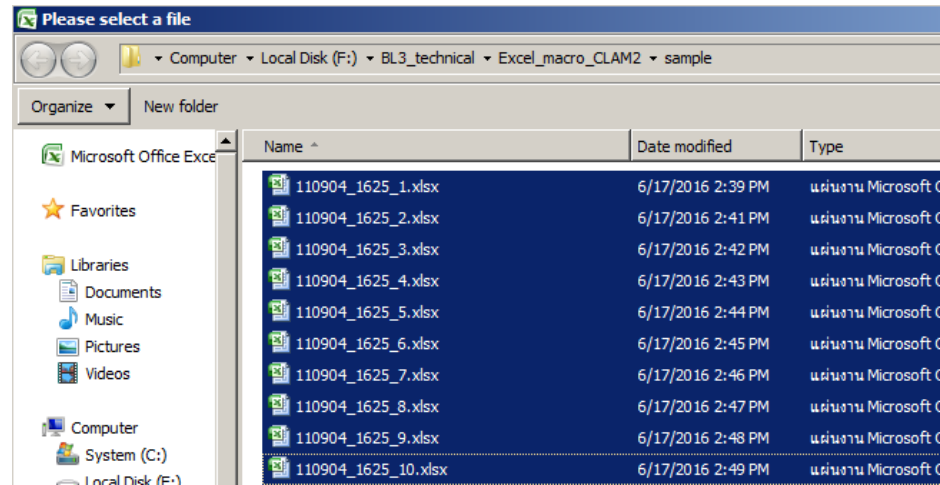

Summary of fitting results

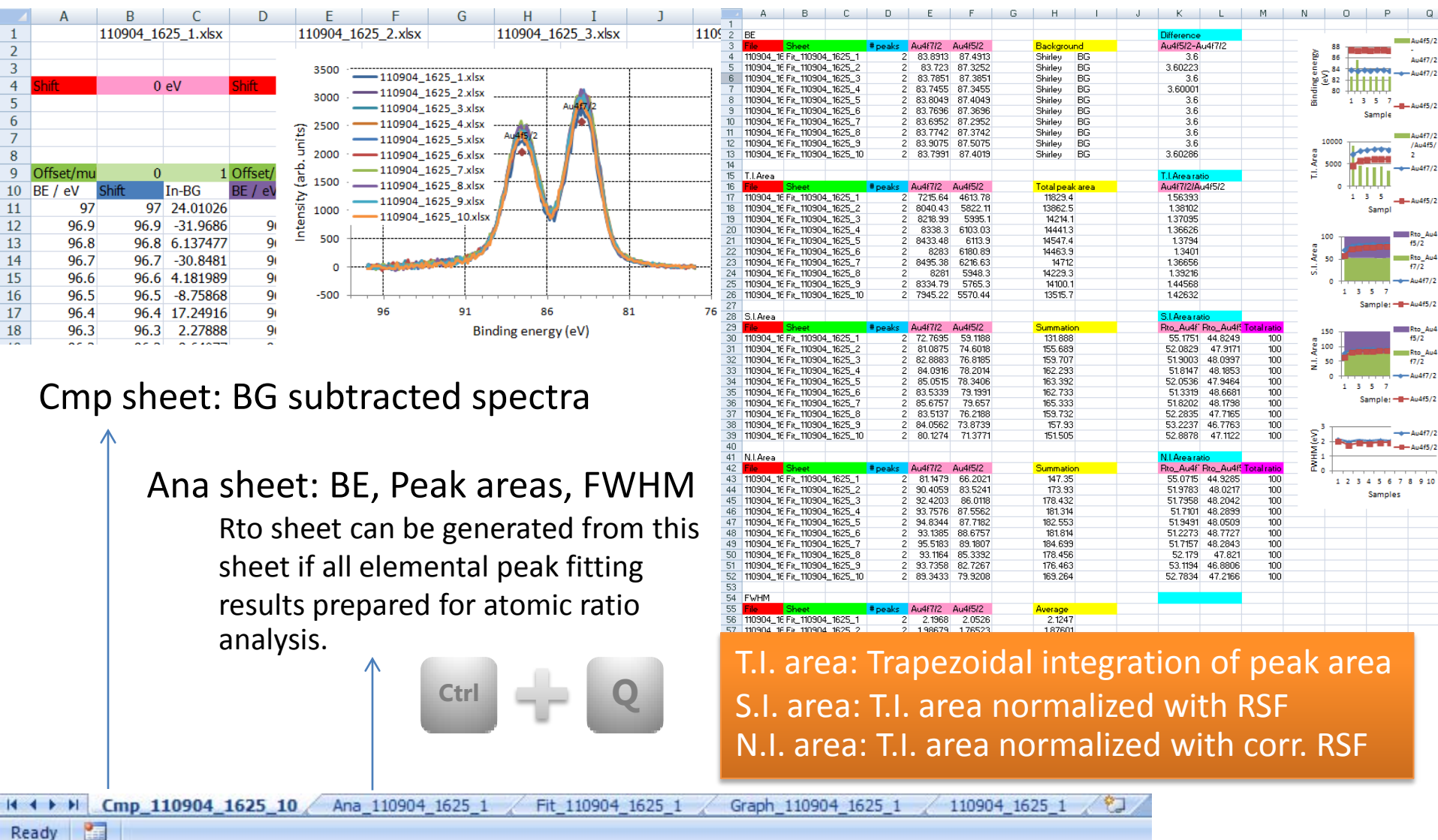

Enjoy!

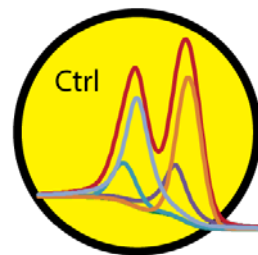

Depth-profile XPS: under review in RSC Adv.

